

# Using VMD: An Introductory Tutorial

UNIT 5.7

Jen Hsin,<sup>1,2</sup> Anton Arkhipov,<sup>1,2</sup> Ying Yin,<sup>1,2</sup> John E. Stone,<sup>2</sup> and Klaus Schulten<sup>1,2</sup>

<sup>1</sup>Department of Physics, University of Illinois at Urbana-Champaign, Urbana, Illinois

<sup>2</sup>Beckman Institute, University of Illinois at Urbana-Champaign, Urbana, Illinois

## ABSTRACT

VMD (Visual Molecular Dynamics) is a molecular visualization and analysis program designed for biological systems such as proteins, nucleic acids, lipid bilayer assemblies, etc. This unit will serve as an introductory VMD tutorial. We will present several step-by-step examples of some of VMD's most popular features, including visualizing molecules in three dimensions with different drawing and coloring methods, rendering publication-quality figures, animating and analyzing the trajectory of a molecular dynamics simulation, scripting in the text-based Tcl/Tk interface, and analyzing both sequence and structure data for proteins. *Curr. Protoc. Bioinform.* 24:5.7.1-5.7.48. © 2008 by John Wiley & Sons, Inc.

Keywords: molecular modeling • molecular dynamics visualization • interactive visualization • animation

## INTRODUCTION

VMD (Visual Molecular Dynamics; Humphrey et al., 1996) is a molecular visualization and analysis program designed for biological systems such as proteins, nucleic acids, lipid bilayer assemblies, etc. It is developed by the Theoretical and Computational Biophysics Group at the University of Illinois at Urbana-Champaign. Among molecular graphics programs, VMD is unique in its ability to efficiently operate on multi-gigabyte molecular dynamics trajectories, its interoperability with a large number of molecular dynamics simulation packages, and its integration of structure and sequence information.


Key features of VMD include methods: (1) general 3D molecular visualization with extensive drawing and coloring methods (e.g., see Fig. 5.7.1); (2) extensive atom selection syntax for choosing subsets of atoms for display; (3) visualization of dynamic molecular data; (4) visualization of volumetric data; (5) support for most molecular data file formats; (6) no limits on the number of atoms, molecules, or trajectory frames, except available memory; (7) molecular analysis commands; (8) rendering high-resolution, publication-quality molecule images; (9) movie making capability; (10) building and preparing systems for molecular dynamics simulations; (11) interactive molecular dynamics simulations; (12) extensions to the Tcl/Python scripting languages; and (13) extensible source code written in C and C++.

This unit will serve as an introductory VMD tutorial. It is impossible to cover all of VMD's capabilities in one unit; instead, we will present several step-by-step examples of VMD's basic features. Topics covered in this tutorial include visualizing molecules in three dimensions with different drawing and coloring methods, rendering publication-quality figures, animating and analyzing the trajectory of a molecular dynamics simulation, scripting in the text-based Tcl/Tk interface, and analyzing both sequence and structure data for proteins.

Modeling  
Structure from  
Sequence

## 5.7.1

Supplement 24

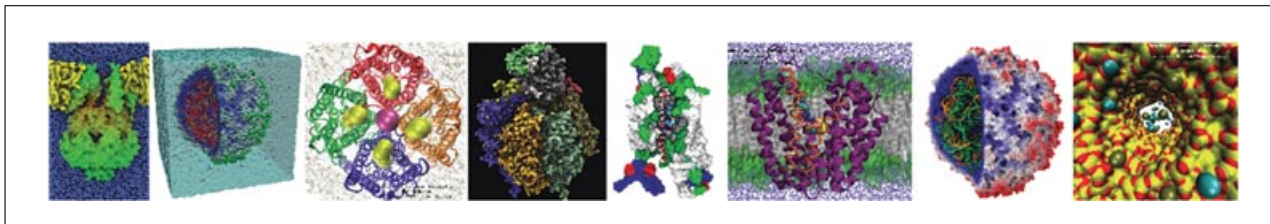

**Figure 5.7.1** Example renderings made with VMD (Cruz-Chu et al., 2006; Freddolino et al., 2006; Yin et al., 2006; Yu et al., 2006; Sotomayor et al., 2007; Wang et al., 2007). For the color version of this figure go to <http://www.currentprotocols.com>.

## DOWNLOADING VMD

Before starting, the current version of VMD needs to be downloaded. This tutorial was written for VMD version 1.8.6. VMD supports all major computer platforms and can be downloaded from the VMD homepage <http://www.ks.uiuc.edu/Research/vmd>. Follow the instructions online to install. Once VMD is installed, to start VMD if using Mac OS X, double-click on the VMD application icon in the Applications directory; if using Linux and SUN, type `vmd` in a terminal window, or if using Windows, select → Start Programs → VMD.

When VMD starts, by default three windows will open: the VMD Main window, the OpenGL Display window, and the VMD Console window (or a Terminal window on a Mac). To end a VMD session, go to the VMD Main window, and choose File → Quit. You can also quit VMD by closing the VMD Console window or the VMD Main window.

## TOPICS AND FILES

This unit contains six sections. Each section acts as an independent tutorial for a specific topic (Working with a Single Molecule, Trajectories and Movie Making, Scripting in VMD, Working with Multiple Molecules, Comparing Protein Structures and Sequences with the MultiSeq Plugin, and Data Analysis in VMD). For readers with no prior experience with VMD, we suggest they work through the sections in the order they are presented. Readers already familiar with the basics of VMD may selectively pursue sections of their interest. Several files have been prepared to accompany this tutorial. You need to download these files at <http://www.currentprotocols.com>.

## WORKING WITH A SINGLE MOLECULE

In this section, the basic functions of VMD will be introduced, starting with loading a molecule, displaying the molecule, and rendering publication-quality molecule images. This section uses the protein ubiquitin as an example molecule. Ubiquitin is a small protein responsible for labeling proteins for degradation, and is found in all eukaryotes with nearly identical sequences and structures.

### *Necessary Resources*

#### *Hardware*

Computer

#### *Software*

VMD, and an image-displaying program

#### *Files*

1ubq.pdb, which can be downloaded at <http://www.currentprotocols.com>

## Loading and Displaying the Molecule

A VMD session usually starts with loading structural information of a molecule into VMD. When VMD loads a molecule, it accesses the information about the names and coordinates of the atoms. Then, one can explore various VMD visualization features to get a nice view of the loaded molecule.

### Loading a molecule

The first step is to load the molecule. The pdb file, 1ubq.pdb (Vijay-Kumar et al., 1987) that contains the atomic coordinates of ubiquitin will be loaded.


1. Start a VMD session. In the VMD Main window, choose File → New Molecule... (Fig. 5.7.2A). The Molecule File Browser window (Fig. 5.7.2B) will appear on the screen.
2. Use the Browse... (Fig. 5.7.2C) button to find the file 1ubq.pdb. When the file is selected, you will be back in the Molecule File Browser window. In order to actually load the file, press Load (Fig 5.7.2D).
3. Now, ubiquitin is shown in the OpenGL Display window. Close the Molecule File Browser window at any time.

*VMD can download a pdb file from the Protein Data Bank (<http://www.pdb.org>) if a network connection is available. Just type the four letter code of the protein in the File Name text entry of the Molecule File Browser window and press the Load button. VMD will download it automatically.*


### Displaying the molecule

In order to see the 3D structure of our protein, the mouse will be used in multiple modes to change the viewpoint. VMD allows users to rotate, scale, and translate the viewpoint of the molecule.


4. In the OpenGL Display, press the left mouse button down and move the mouse. Explore what happens. This is the rotation mode of the mouse and allows for rotation of the molecule around an axis parallel to the screen (Fig. 5.7.3A).

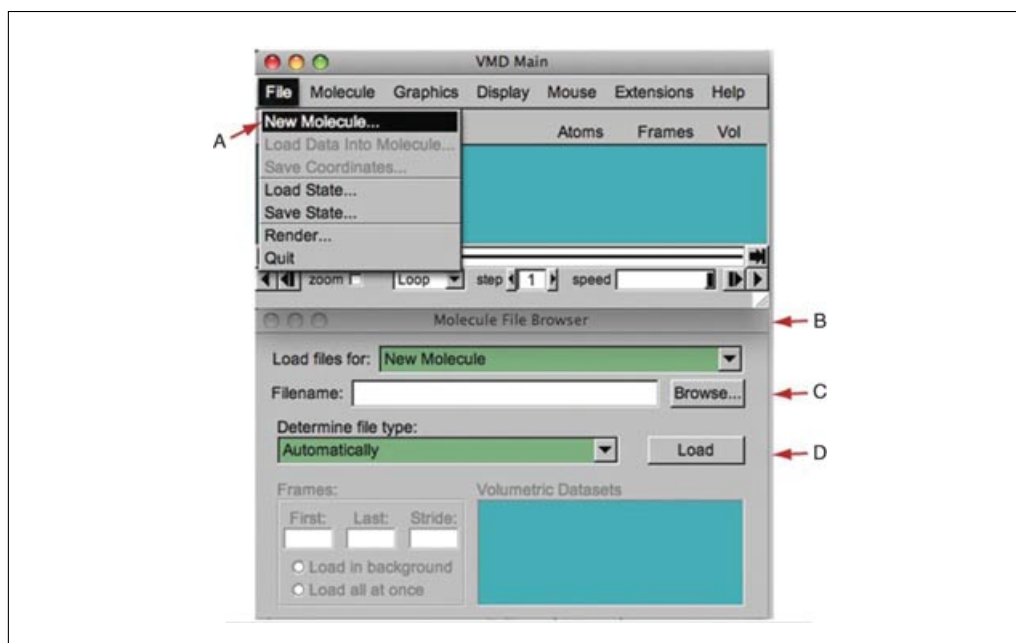

**Figure 5.7.2** Loading a molecule.

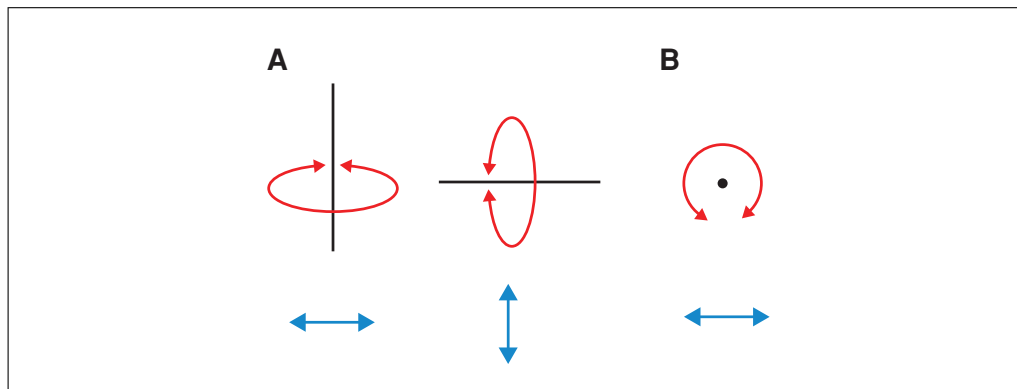

**Figure 5.7.3** Rotational modes. (A) Rotation axes when holding down the left mouse key. (B) The rotation axes when holding down the right mouse key. For the color version of this figure go to <http://www.currentprotocols.com>.

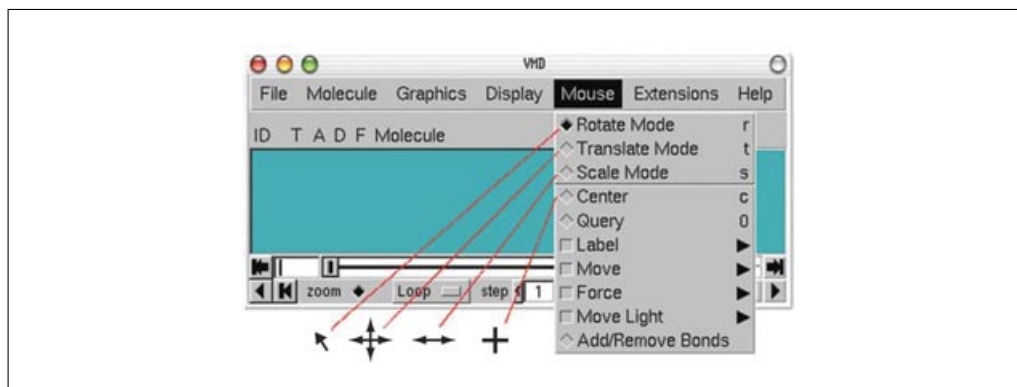

**Figure 5.7.4** Mouse modes and their characteristic cursors.

5. Holding down the right mouse button and repeating the previous step will cause rotation around an axis perpendicular to the screen (Fig. 5.7.3B).

*For Mac users who have a single-button mouse or a trackpad, the right mouse button is equivalent to holding down the command key while pressing the mouse/trackpad button.*


6. In the VMD Main window, look at the Mouse menu (Fig. 5.7.4). Here, the user is able to switch the mouse mode from Rotation to Translation or Scale modes.
7. Choose the Translation mode and go back to the OpenGL Display. It is now possible to move the molecule around when you hold the left mouse button down.
8. Go back to the Mouse menu and choose the Scale mode this time. This will allow the user to zoom in or out by moving the mouse horizontally while holding down the left mouse button.

*It should be noted that these actions performed with the mouse only change the viewpoint and do not change the actual coordinates of the molecule's atoms.*

*Also note that each mouse mode has its own characteristic cursor and its own shortcut key (r: Rotate, t: Translate, s: Scale). When the OpenGL Display window is the active window, these shortcut keys can be used instead of the Mouse menu to change the mouse mode.*

*Another useful option is the Mouse → Center menu item. It allows you to specify the point around which rotations are done.*

9. Select the Center menu item and pick one atom at one of the ends of the protein; the cursor should display a cross.

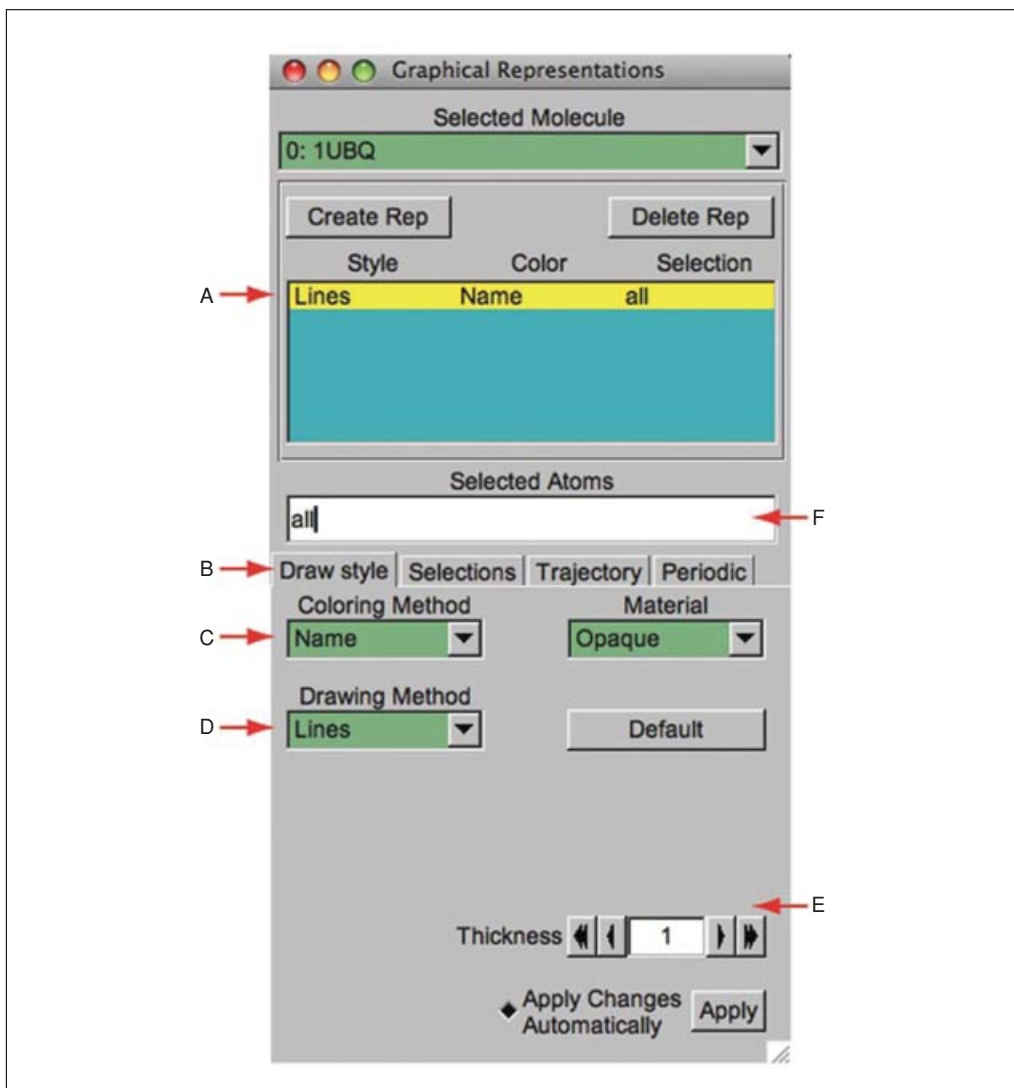


**Figure 5.7.5** The Graphical Representations window. (A) List of representations, (B) the tabs for Draw Style, Selections, Trajectory, and Periodic, (C) Coloring Method pull-down menu, (D) Drawing Method pull-down menu, (E) user-adjustable parameters for different drawing methods, and (F) selection text entry box.

10. Now, press  $r$ , and rotate the molecule with the mouse and see how the molecule moves around the selected point.
11. In the VMD Main window, select the Display → Reset View menu item to return to the default view. You can also reset the view by pressing the “=” key when you are in the OpenGL Display window.

### **Graphical representations**

VMD can display molecules in various ways by setting the Graphical Representations window shown in Figure 5.7.5. Each representation is defined by four main parameters: the selection of atoms included in the representation, the drawing style, the coloring method, and the material. The selection determines which part of the molecule is drawn; the drawing method defines which graphical representation is used; the coloring method gives the color of each part of the representation; and the material determines the effects of lighting, shading, and transparency on the representation. Let us first explore different drawing styles.

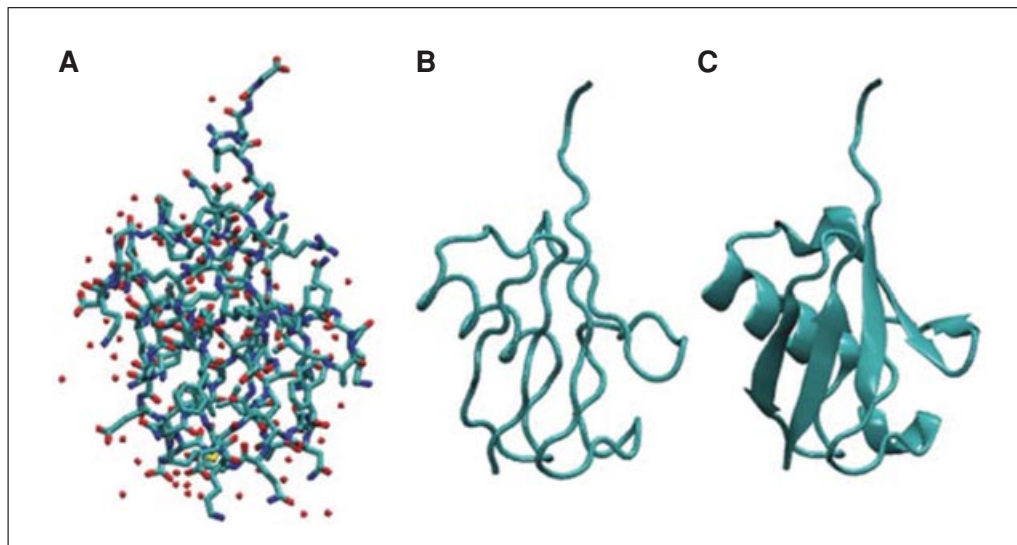

**Figure 5.7.6** (A) Licorice, (B) Tube, and (C) NewCartoon representations of ubiquitin. For the color version of this figure go to <http://www.currentprotocols.com>.

#### Exploring different drawing styles

12. In the VMD Main window, choose the Graphics → Representations... menu item. A window called Graphical Representations will appear and the current default representation will be highlighted in yellow (Fig. 5.7.5A).
13. In the Draw Style tab (Fig. 5.7.5B), change the style (Fig. 5.7.5D) and color (Fig. 5.7.5C) of the representation. Here, we will focus on the drawing style (the default is Lines).
14. Each Drawing Method has its own parameters. For instance, change the thickness of the lines by using the controls on the lower right-hand-side corner (Fig. 5.7.5E) of the Graphical Representations window.
15. Click on the Drawing Method (Fig. 5.7.5D) to see a list of options. Choose VDW (van der Waals); each atom is now represented by a sphere scaled to its van der Waals radius, allowing the user to see the volumetric distribution of the protein.
16. When choosing VDW as the drawing method, two new controls will show up in the lower right-hand-side corner. Use these controls to change the Sphere Scale to 0.5 and the Sphere Resolution to 13. Note that the higher the resolution, the slower the display of the molecule will be.
17. Press the Default button. This returns the screen to the default properties of the chosen drawing method.

*Other popular representations include CPK and Licorice. In CPK, like in old chemistry ball and stick kits, each atom is represented by a sphere and each bond is represented by a thin cylinder (radius and resolution of both the sphere and the cylinder can be modified). The Licorice drawing method also represents each atom as a sphere and each bond as a cylinder, but the sphere and the cylinder have the same radii.*

#### Using the Tube style drawing method

The previous representations visualize micromolecular details of the protein by displaying every single atom. More general structural properties can be demonstrated better by using more abstract drawing methods.

18. Choose the Tube style under Drawing Method, which shows the backbone of the protein. Set the Radius to 0.8. The result should be similar to Figure 5.7.6.

### *Using the NewCartoon drawing method*

The last drawing method described here is NewCartoon. It gives a simplified representation of a protein based on its secondary structure. Helices are drawn as coiled ribbons,  $\beta$ -sheets as solid, flat arrows, and all other structures as a tube. This is probably the most popular drawing method to view the overall architecture of a protein.

19. In the Graphical Representations window, choose Drawing Method  $\rightarrow$  NewCartoon. The helices,  $\beta$ -sheets, and coils of the protein can now be easily identified.

*Ubiquitin has three and one half turns of  $\alpha$ -helix (residues 23 to 34, three of them hydrophobic), one short piece of  $3_{10}$ -helix (residues 56 to 59) and a mixed  $\beta$ -sheet with five strands (residues 1 to 7, 10 to 17, 40 to 45, 48 to 50, and 64 to 72), and seven reverse turns. VMD uses the program STRIDE (Frishman and Argos, 1995) to compute the secondary structure according to a heuristic algorithm.*

### *Exploring different coloring methods*

In this series of steps, different coloring methods are explored.

20. In the Graphical Representations window, the default coloring method is Coloring Method  $\rightarrow$  Name. In this coloring method, choose a drawing method that shows individual atoms: each atom will have a different color, i.e., O is red, N is blue, C is cyan, and S is yellow.
21. Choose Coloring Method  $\rightarrow$  ResType (Fig. 5.7.5C). This allows nonpolar residues (white) to be distinguished from basic residues (blue), acidic residues (red), and polar residues (green).
22. Select Coloring Method  $\rightarrow$  Structure (Fig. 5.7.5C) and confirm that the NewCartoon representation displays colors consistent with secondary structure.

### *Displaying different selections*


To display only parts of the molecule of interest, one can specify their selection in the Graphical Representations window (Fig. 5.7.5F).

23. In the Graphical Representations window, there is a Selected Atoms text entry (Fig. 5.7.5F). Delete the word `all`, type `helix`, and press the Apply button or hit the Enter/return key (remember to do this whenever a selection is changed). VMD will show just the helices present in the molecule.
24. In the Graphical Representations window, choose the Selections tab (Fig. 5.7.7A). In the section Singlewords (Fig. 5.7.7B), a list of possible selections that can be entered is provided.

*Combinations of Boolean operators can also be used when writing a selection.*

25. In order to see the molecule without helices and  $\beta$ -sheets, type the following in the Selected Atoms field: `(not helix) and (not betasheet)`. Remember to press the Apply button or hit the Enter/return key.
26. In the section Keyword (Fig. 5.7.7C) of the Selections tab, the properties that can be used to select parts of a molecule are listed along with their possible values. Look at possible values of the keyword "resname" (Fig. 5.7.7D).
27. Display all the lysines and glycines present in the protein by typing `(resname LYS) or (resname GLY)` in the Selected Atoms field.

*Lysines play a fundamental role in the configuration of polyubiquitin chains.*

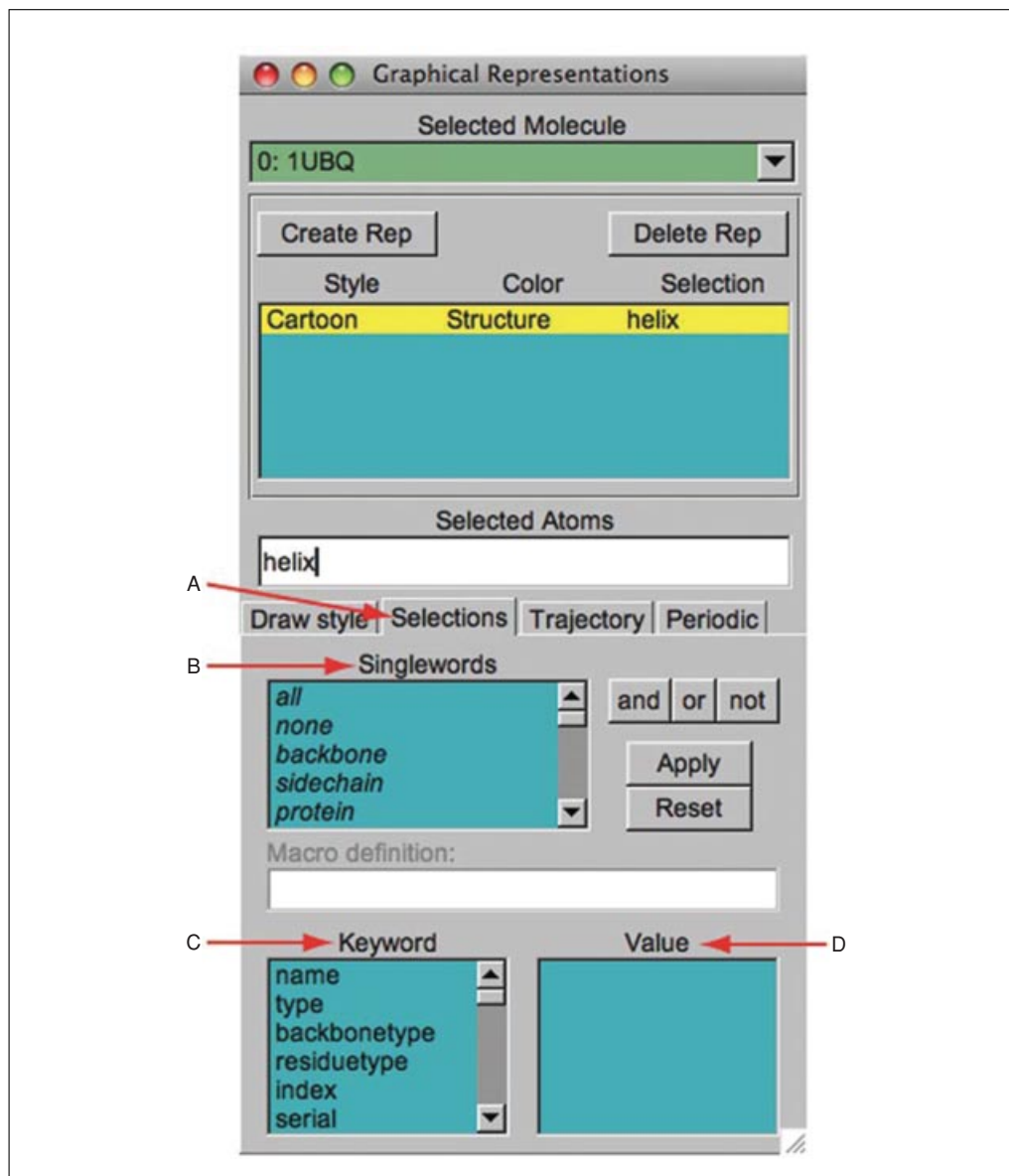

**Figure 5.7.7** Graphical Representations window and the (A) Selections tab, (B) list of Singlewords, (C) list of Keywords, and (D) Value box that displays possible choices for a given keyword.


28. Change the current representation's Drawing Method to CPK and the Coloring Method to ResName in the Draw style tab. In the screen, the different lysines and glycines will be visible.
29. In the Selected Atoms text field, entry type `water`. Choose Coloring Method → Name. The 58 water molecules present in the system now appear (in fact only their oxygen atoms).
30. In order to see which water molecules are closer to the protein, use the command `within`. Type `water` and `within 3` of `protein` for Selected Atoms in the text field.

*This selects all the water molecules that are within a distance of 3 Å of the protein.*

31. Finally, try typing in the Selected Atoms field the selections shown in the first column of Table 5.7.1. Each of these selections will show the protein or part of the protein as explained in the second column of Table 5.7.1.

**Table 5.7.1** Examples of Atom Selections

Selection	Action
Protein	Shows the protein
resid 1	The first residue
(resid 1 76) and (not water)	The first and last residues
(resid 23 to 34) and (protein)	The $\alpha$ -helix

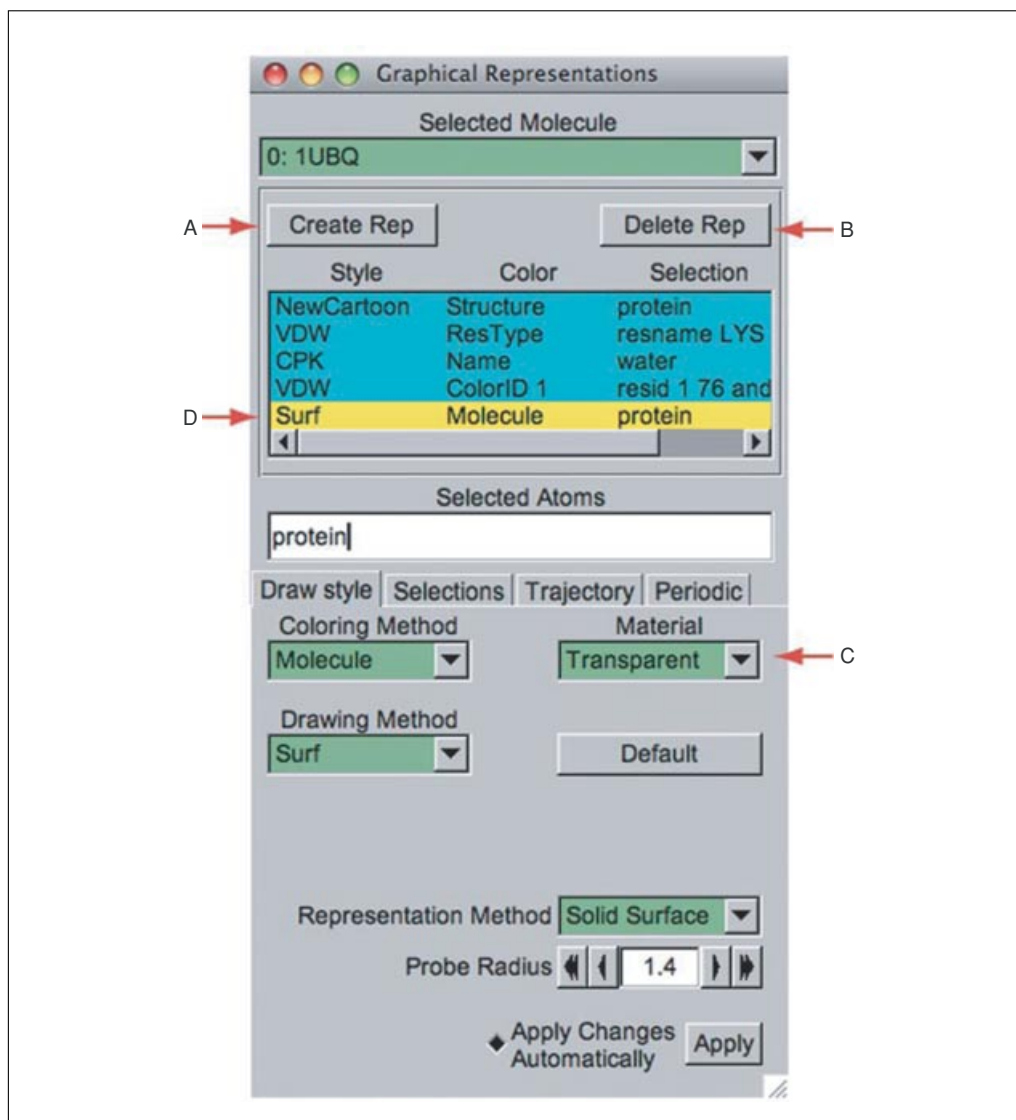

**Figure 5.7.8** Multiple Representations of ubiquitin. Representations can be either created or deleted using the (A) Create Rep and (B) Delete Rep buttons. Screen also shows (C) the Material pull-down menu and (D) list of representations.

### *Creating multiple representations*

The button Create Rep (Fig. 5.7.8A) in the Graphical Representations window allows creation of multiple representations. Therefore, users can have a mixture of different selections with different styles and colors, all displayed at the same time.

- For the current representation, in the Selected Atoms field type `protein`, set the Drawing Method to NewCartoon and the Coloring Method to Structure in the Draw style tab.

**Table 5.7.2** Examples of Representations

Selection	Coloring method	Drawing method
Water	Name	CPK
resid 1 76 and name CA	ColorID 1	VDW

33. Press the Create Rep button (Fig. 5.7.8A). A new representation will be created.
34. Modify the new representation to get VDW as the Drawing Method, ResType as the Coloring Method, and `resname LYS` as the current selection.
35. Repeating the previous procedure, create the following two new representations in Table 5.7.2. These two representations show water molecules and the C $\alpha$  atoms of the first and last residues of the protein.
36. Create the last representation by pressing the Create Rep button again. Select Drawing Method  $\rightarrow$  Surf for drawing method, Coloring Method  $\rightarrow$  Molecule for coloring method, and type `protein` in the Selected Atoms field. For this last representation, choose Transparent in the Material pull-down menu (Fig. 5.7.8C). This representation shows the protein's volumetric surface in transparent.


*Note that you can select and modify different representations you have created by clicking on a representation to highlight it in yellow. Also, each representation can be switched on/off by double-clicking on it. To delete a representation, highlight it and then click on the Delete Rep button (Fig. 5.7.8B). At the end of this section, the Graphical Representations window should look like Figure 5.7.8.*

### Sequence viewer extension

When dealing with a protein for the first time, it is very useful to be able to find and display different amino acids quickly. The sequence viewer extension allows viewing of the protein sequence, as well as to easily pick and display one or more residues of interest.

37. In the VMD Main window, choose the Extension  $\rightarrow$  Analysis  $\rightarrow$  Sequence Viewer menu item. A window (Fig. 5.7.9A) with a list of the amino acids (Fig. 5.7.9E) and their properties (Figs. 5.7.9B through 5.7.9C) will appear on the screen.
38. With the mouse, try clicking on different residues in the list (Fig. 5.7.9E) and see how they are highlighted. In addition, the highlighted residue will appear in the OpenGL Display window in yellow and rendered in the bond drawing method, so its location within the protein can be visualized easily.
39. Use the Zoom controls (Fig. 5.7.9F) to display the entire list of residues in the window. This is particularly useful for larger proteins.
40. Pick multiple residues by holding the shift key and clicking on the mouse button (Fig. 5.7.9E).
41. Look at the Graphical Representations window; a new representation with the residues that have been selected using the Sequence Viewer Extension should be shown. Modify, hide, or delete this representation similar to the steps described above.

*Information about residues is color-coded (Fig. 5.7.9D) in columns and obtained from STRIDE. The B-value column (Fig. 5.7.9B) shows the B-value field (temperature factor) often provided in `pdb` files. The "struct" column shows secondary structure (Fig. 5.7.9D), where each letter corresponds to a secondary structure, listed in Table 5.7.3.*

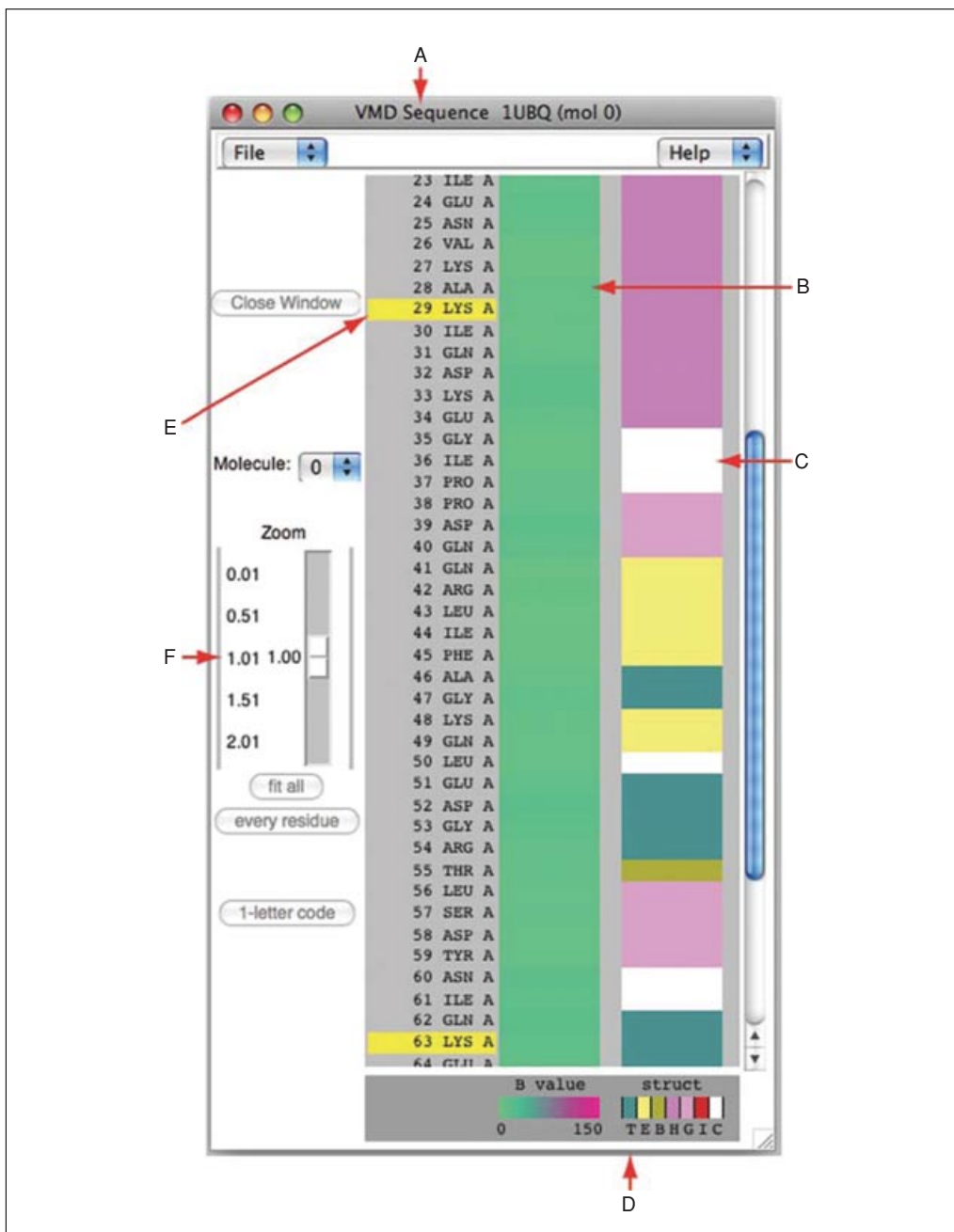

**Figure 5.7.9** (A) The VMD Sequence window displays properties of the protein sequence, including (B) the B-value and (C) the secondary structure, denoted by (D) the color codes. (E) The list of residues is displayed, with the selected residues highlighted in yellow. (F) Zoom controls are also shown in the window. For the color version of this figure go to <http://www.currentprotocols.com>.

### ***Saving your work***

The viewpoints and representations created using VMD can be saved as a VMD state. This VMD state contains all the information needed to reproduce the same VMD session.

42. Go to the OpenGL Display window; use the mouse to find a nice view of the protein. We will save this viewpoint using the VMD ViewMaster.
43. In the VMD Main window (Fig. 5.7.2), select Extension → Visualization → ViewMaster. This will open the VMD ViewMaster window.
44. In the VMD ViewMaster window, click on the Create New button. The OpenGL Display viewpoint has now been saved.

**Table 5.7.3** Secondary Structure Codes Used by STRIDE

Letter code	Secondary structure
T	Turn
E	Extended conformation ( $\beta$ -sheets)
B	Isolated bridge
H	Alpha helix
G	3-10 helix
I	Pi helix
C	Coil

45. Go back to the OpenGL Display window and use the mouse to find another nice view. If desired, you can add/delete/modify a representation in the Graphical Representations window. When a good view has been found, save it by returning to the VMD ViewMaster window and clicking on the Create New button.
46. Create as many views as desired by repeating the previous step. All of the viewpoints are displayed as thumbnails in the VMD ViewMaster window. A previously saved viewpoint can be opened by clicking on its thumbnail.
47. To save the entire VMD session, in the VMD Main window, choose the File  $\rightarrow$  Save State menu item. Type an appropriate name (e.g., `myfirststate.vmd`) and save it.

*The VMD state file `myfirststate.vmd` contains all the information needed to restore a VMD session, including the viewpoints and the representations.*

*To load a saved VMD state, start a new VMD session and in the VMD Main window choose File  $\rightarrow$  Load State.*

48. Quit VMD.

## BASIC PROTOCOL 2

### The Basics of VMD Figure Rendering

One of VMD's many strengths is its ability to render high-resolution, publication-quality molecule images. In this section, we will introduce some basic concepts of figure rendering in VMD.

#### *Setting the display background*

Before rendering a figure, make sure that the OpenGL Display background is set up the way you want. Nearly all aspects of the OpenGL Display are user-adjustable, including the background color.

1. Start a new VMD session (Basic Protocol 1) and load the `1ubq.pdb` file.
2. In the VMD Main window, choose Graphics  $\rightarrow$  Colors... The Color Controls window should show up. Look through the Categories list. All display colors, for example, the colors of different atoms when colored by name, are set here.
3. In Categories, select Display. In Names, select Background. Finally, choose "8 white" in Colors. The OpenGL Display should now have a white background.
4. When making a figure, we often do not want to include the axes. To turn off the axes, select Display  $\rightarrow$  Axes  $\rightarrow$  Off in the VMD Main window.

### ***Increasing geometric resolution***

All VMD objects are drawn with an adjustable resolution, allowing users to balance fineness of detail with drawing speed.

5. Open the Graphical Representation window via Graphics → Representations... in the VMD Main menu. Modify the default representation to show just the protein, and display it using the VDW drawing method.
6. Zoom in on one or two of the atoms by using Mouse Scale → Mode (shortcut s).


*You might notice that as you zoom into an atom closer and closer, the atom might be cut off by an invisible clipping plane, which makes it difficult to focus on just one atom. This is an OpenGL feature. You can move the clipping plane closer to you by doing the following: switch your mouse mode to the Translate mode, either by pressing the shortcut key “t” in the OpenGL window or by selecting Mouse → Translate Mode, and dragging your mouse in the OpenGL window while holding down the right mouse key. You can now move the clipping plane closer to you, or away from you. If this does not work, here is an alternative way: in the VMD Main window, choose Display → Display Settings...; in the Display Settings window that shows up you can see that many OpenGL options are adjustable; decrease the value for Near Clip, which will move the OpenGL clipping closer, allowing you to zoom in on individual atoms without clipping them off.*

7. Notice that with the default resolution setting, the “spherical” atoms are not looking very spherical. In the Graphical Representations window, click on the representation you set up before for the protein to highlight it in yellow. Try adjusting the Sphere Resolution setting to something higher, and see what a difference it can make (Fig. 5.7.10).

*Most of the drawing methods have a geometric resolution setting. Try a few different drawing methods and see how their resolutions can be easily increased. When producing images, the resolution can be raised until it stops making a visible difference.*

### ***Colors and materials***

8. There is a Material menu in the Graphical Representations window (which by default is set to Opaque material). Choose the protein representation you made before, and experiment with the different materials in the Material menu.

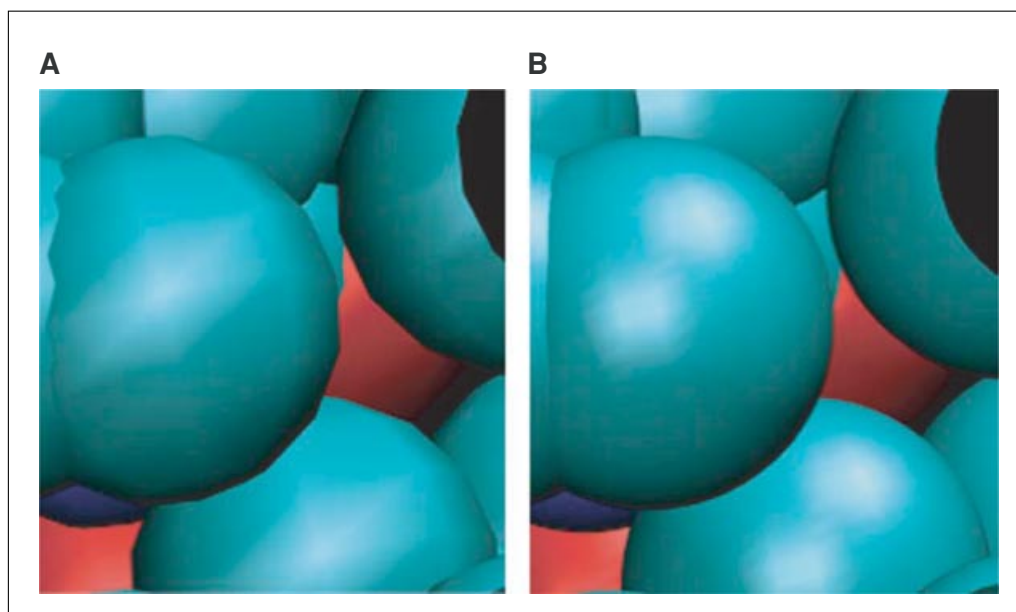

**Figure 5.7.10** The effect of the resolution setting. (A) Low resolution: Sphere Resolution set to 8. (B) High resolution: Sphere Resolution set to 28.


9. Besides the predefined materials in the Material menu, VMD also allows users to create their own materials. To make a new material, in the VMD Main window choose Graphics → Materials. . .
10. In the Materials window that appears, you will see a list of the materials you just tried out, and their adjustable settings. Click the Create New button. A new material, Material 12, will be created. Give it the settings listed in Table 5.7.4.
11. Go back to the Graphical Representations window. In the Material menu, Material 12 is now on the list. Try using Material 12 for a representation and see what it looks like. You can also rename the materials in the Material menu.

*Now is a good time to try out the GLSL Render Mode, if your computer supports it. In the VMD Main window, choose Display → Rendermode → GLSL. This mode uses your 3D graphics card to render the scene with real-time ray-tracing of spheres and alpha-blended transparency, and can improve the visualization of transparent materials. See Figure 5.7.11 for example renderings made in GLSL mode.*

12. If your computer supports GLSL Render Mode, you can try to reproduce Figure 5.7.11. First, turn on the GLSL rendering mode by selecting Display → Rendermode → GLSL in the VMD Main window.
13. Modify Material 12 to be more transparent by entering the values listed in Table 5.7.5 in the Materials window.

**Table 5.7.4** Example of a User-Defined Material

Setting	Value
Ambient	0.30
Diffuse	0.30
Specular	0.90
Shininess	0.50
Opacity	0.95

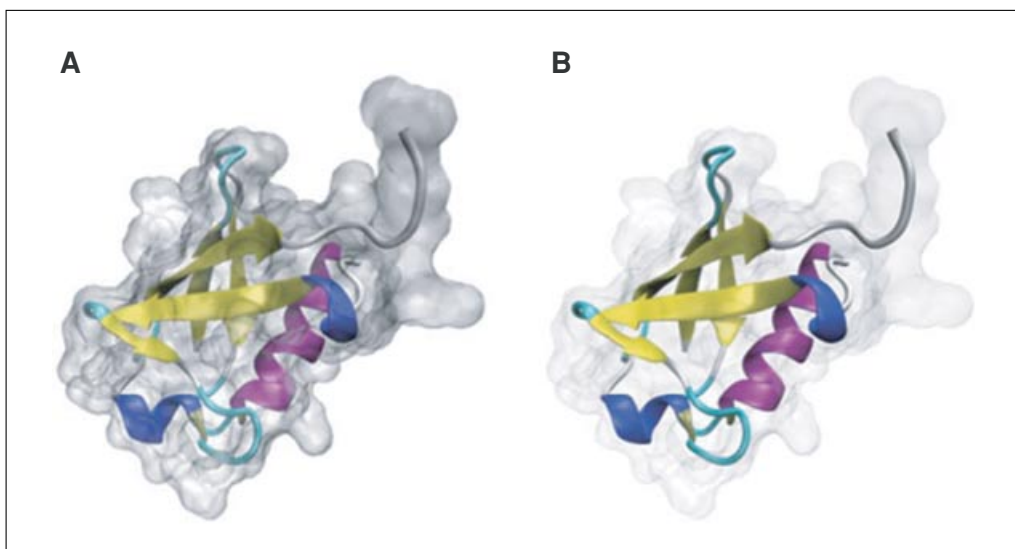

**Figure 5.7.11** Examples of different material settings. **(A)** The default transparent material, rendered in GLSL mode. **(B)** A user-defined material with high transparency, also rendered in GLSG mode. For the color version of this figure go to <http://www.currentprotocols.com>.

**Table 5.7.5** Example of a More Transparent Material

Setting	Value
Ambient	0.30
Diffuse	0.50
Specular	0.87
Shininess	0.85
Opacity	0.11

**Table 5.7.6** Example of Representations Drawn with Different Materials

Selection	Coloring method	Drawing method	Material
protein	Structure	NewCartoon	Opaque
protein	ColorID→8 white	Surf	Material 12

- Hide all of the current representations and create the two representations listed in Table 5.7.6.

### ***Depth perception***

Since the molecular systems are three-dimensional, VMD has multiple ways of representing the third dimension. In this section, how to use VMD to enhance or hide depth perception is discussed.

- The first thing to consider is the projection mode. In the VMD Main window, click the Display menu. Here, we can choose either Perspective or Orthographic in the drop-down menu. Try switching between Perspective or Orthographic projection modes and see the difference (Fig. 5.7.12).


*In perspective mode, things closer to the camera appear larger. Perspective projection provides strong size-based visual depth cues, but the displayed image will not preserve scale relationships or parallelism of lines, and objects very close to the camera may appear distorted. Orthographic projection preserves scale and parallelism relationships between objects in the displayed image, but greatly reduces depth perception. Hence, orthographic mode tends to be more useful for analysis, because alignment is easy to see, while perspective mode is often used for producing figures and stereo images.*

*Another way VMD can represent depth is through so-called “depth cueing.” Depth cueing is used to enhance three-dimensional perception of molecular structures, particularly with orthographic projections.*


- Choose Display → Depth Cueing in the VMD Main window.

*When depth cueing is enabled, objects further from the camera are blended into the background. Depth cueing settings are found in Display → Display Settings. . . . Here one can choose the functional dependence of the shading on distance, as well as some parameters for this function. To see the depth cueing effect better, you might want to hide the representation with the Surf drawing method.*

- Finally, VMD can also produce stereo images. In the VMD Main window, look at the Display → Stereo menu, showing many different choices. Choose SideBySide (remember to return to Perspective mode for a better result). The result should look like Figure 5.7.13.
- Turn off stereo image by selecting Display → Stereo → Off in the VMD Main window. Also, turn off depth cueing by unselecting the Display → Depth Cueing checkbox in the VMD Main window.

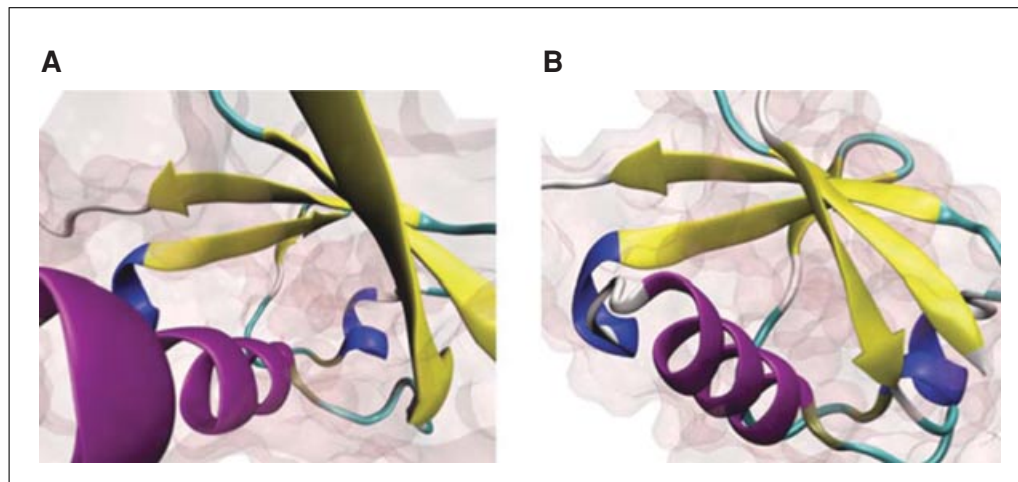

**Figure 5.7.12** Comparison of the (A) perspective and (B) orthographic projection modes. For the color version of this figure go to <http://www.currentprotocols.com>.

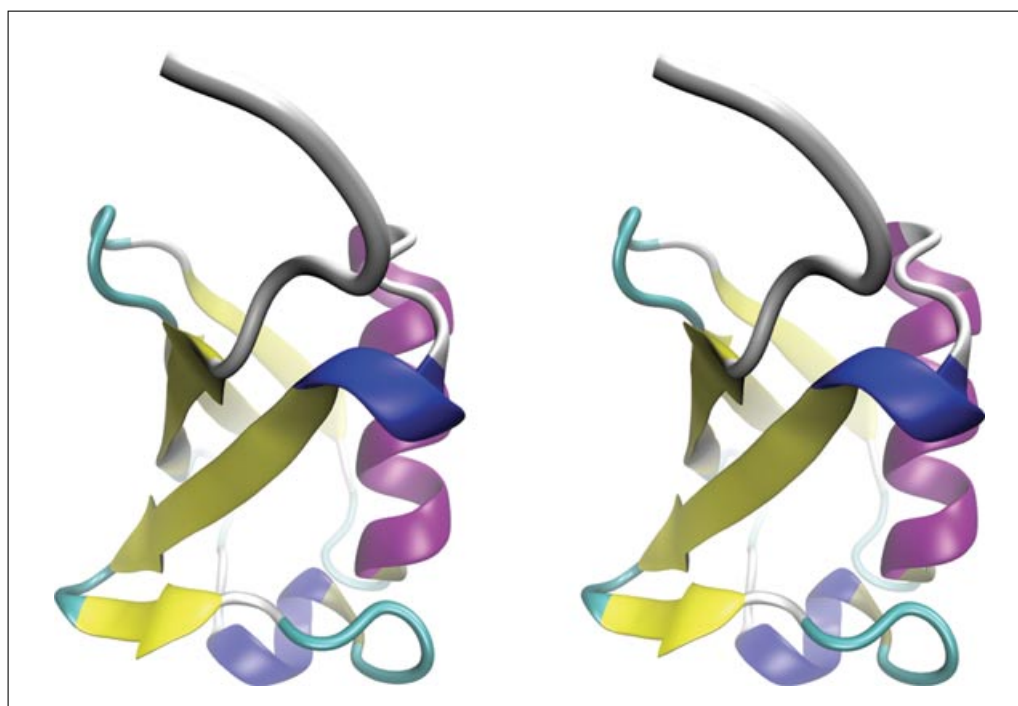

**Figure 5.7.13** Stereo image of the ubiquitin protein. Shown here with Cue Mode = Linear, Cue Start = 1.5, and Cue End = 2.75. To view the stereo image, use the “wall-eyed” method: hold the page close to eyes, and shift the focus beyond the page until the two images overlap to form a three-dimensional object. If this is difficult, try scaling down the figure to a smaller size. This will make viewing easier. For the color version of this figure go to <http://www.currentprotocols.com>.

### **Rendering**

By now, we have seen some techniques for producing nice views and representations of the molecule loaded in VMD. Now, we will explore the use of the VMD built-in snapshot feature and external rendering programs to produce high-quality images of your molecule. The “snapshot” renderer saves the on-screen image in the OpenGL window and is adequate for use in presentations, movies, and small figures. When one desires higher-quality images, renderers such as Tachyon and POV-Ray are better choices.

19. Hide or delete all previous representations, and create the four new representations listed in Table 5.7.7.

**Table 5.7.7** Example Representations

Selection	Coloring method	Drawing style	Material
protein and not resid 72 to 76	Structure	NewCartoon	Opaque
protein and helix and name CA	ColorID→8	Surf	Material 12
resname GLY and not resid 72 to 76	ColorID→7	VDW	Opaque
resname LYS	ColorID→18	Licorice	Opaque


20. Once you have the scene set the way you like it in the OpenGL window, simply choose File → Render. . . in the VMD Main window. The File Render Controls window will appear on the screen.

21. The File Render Controls allows you to choose which renderer you want to use and the file name for your image. Select “snapshot” for the rendering method, type in a filename of your choice, and click Start Rendering.

22. If you are using a Mac or a Linux machine, an image-processing application might open automatically that shows you the molecule you have just rendered using Snapshot. If this is not the case, use any image-processing application to take a look at the image file. Close the application when you are done to continue using VMD.

*The snapshot renderer saves exactly what is showing in your OpenGL display window—in fact, if another window overlaps the display window, it may distort the overlapped region of the image.*

23. Try to render again using different rendering methods, particularly TachyonInternal and POV3 (see Fig. 5.7.14 for an example POV3 rendering). Compare the quality of the images created by different renderers.

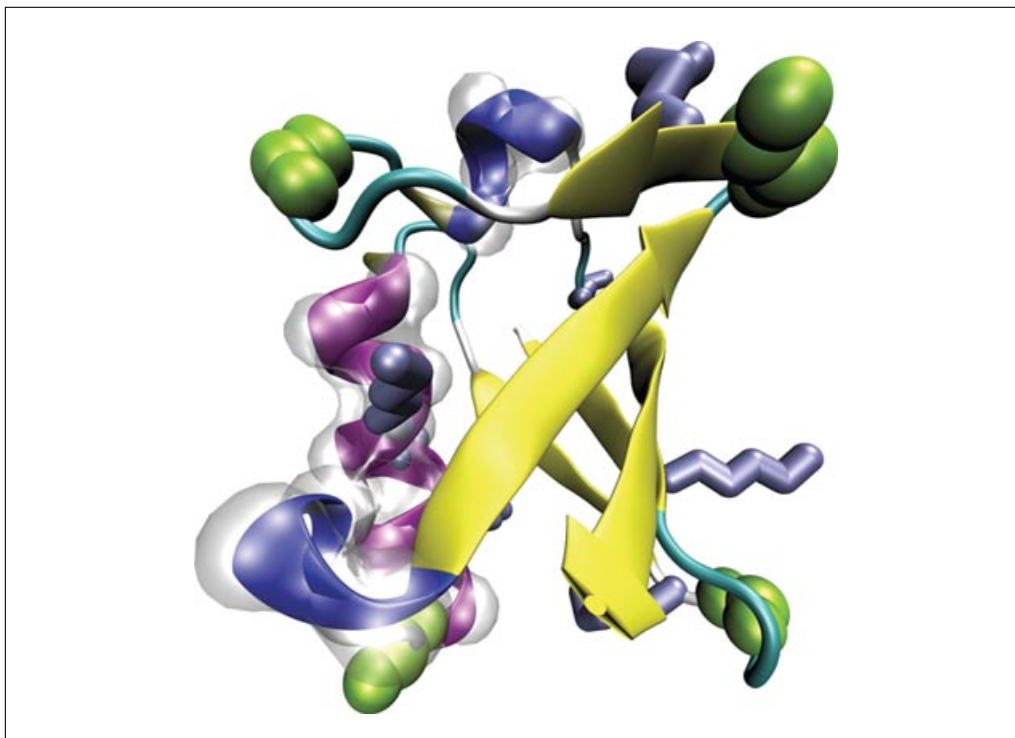

**Figure 5.7.14** Example of a POV3 rendering. For the color version of this figure go to <http://www.currentprotocols.com>.

*The other renderers (e.g., POV3 and Tachyon) reprocess everything, so it may not look exactly as it does in the OpenGL window. In particular, they do not “clip,” or hide, objects very near the camera. If you select Display → Display Settings... in the VMD Main window, you can set Near Clip to 0.01 to get a better idea of what will appear in your rendering.*

24. Quit VMD.

## WORKING WITH TRAJECTORIES AND MAKING MOVIES

Time-evolving coordinates of a system are called trajectories. They are most commonly obtained from simulations of molecular systems, but can also be generated by other means and for different purposes. Upon loading a trajectory into VMD, one can see a movie of how the system evolves in time and analyze various features throughout the trajectory. This section will introduce the basics of working with trajectory data in VMD. You will also learn how to analyze trajectory data in Basic Protocols 14, 15, and 16.

### *Necessary Resources*

#### *Hardware*

Computer

#### *Software*

VMD, and a movie player program

#### *Files*

`ubiquitin.psf` and `pulling.dcd`, which can be downloaded from <http://www.currentprotocols.com>

### Working with Trajectories

Trajectory files are commonly binary files that contain several sets of coordinates for the system. Each set of coordinates corresponds to one frame in time. An example of a trajectory file is a DCD file generated by the molecular dynamics program NAMD (Phillips et al., 2005).


#### *Load trajectories*

Trajectory files do not contain information of the system contained in the protein structure files (PSF). Therefore, we first need to load the PSF file, and then add the trajectory data to this file.

1. Start a new VMD session. In the VMD Main window, select File → New Molecule. . . The Molecule File Browser window will appear on your screen.
2. Use the Browse. . . button to find the file `ubiquitin.psf`. When you select this file, you will be back in the Molecule File Browser window. Press the Load button to load the molecule.
3. In the Molecule File Browser window, make sure that `ubiquitin.psf` is selected in the “Load files for:” pull-down menu on top, and click on the Browse button. Browse for `pulling.dcd`.

*Note the options available in the Molecule File Browser window: one can load trajectories starting and finishing at chosen frames, and adjust the stride between the loaded frames. Leave the default settings so that the whole trajectory is loaded.*

4. Click on the Load button in the Molecule File Browser window.

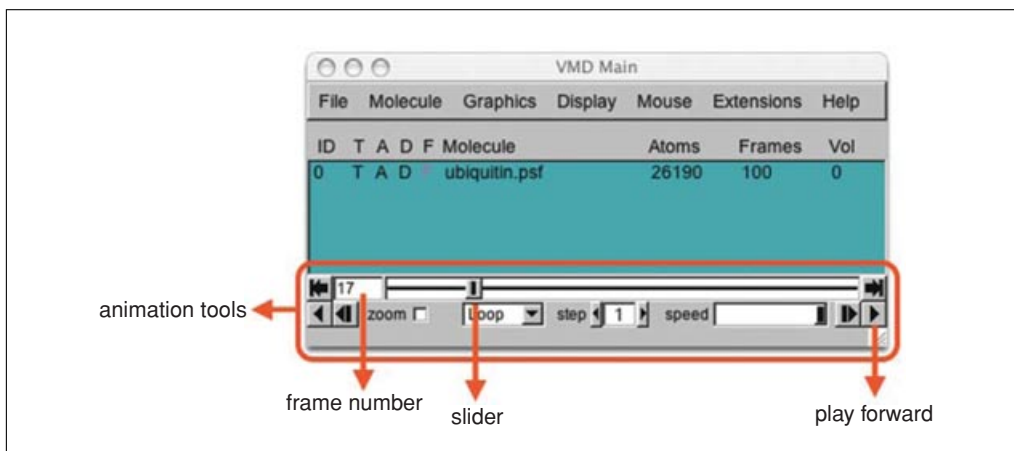

**Figure 5.7.15** Animation tools in the VMD main menu. The tools allow one to go over frames of the trajectory (e.g., using the “slider”) and to play a movie of the trajectory in various modes (Once, Loop, or Rock) and at an adjustable speed.

*You will be able to see the frames as they are loaded into the molecule in the OpenGL window. After the trajectory finishes loading, you will be looking at the last frame of your trajectory. To go to the beginning, use the animation tools at the lower part of the VMD Main menu (see Fig. 5.7.15).*

5. Close the Molecule File Browser window.
6. For a convenient visualization of the protein, choose Graphics → Representations in the VMD Main menu. In the Selected Atoms text field, type `protein` and hit Enter on your keyboard; in the Drawing Method, select NewCartoon; in the Coloring Method, select Structure.

*The trajectory you just loaded is a simulation of an AFM (Atomic Force Microscopy) experiment pulling on a single ubiquitin molecule, performed using the Steered Molecular Dynamics (SMD) method (Isralewitz et al., 2001). We are looking at the behavior of the protein as it unfolds while being pulled from one end, with the other end constrained to its original position. Each frame corresponds to 10 picoseconds of simulation time. Ubiquitin has many functions in the cell, and it is currently believed that some of these functions depend on the protein’s elastic properties, which can be probed in AFM pulling experiments. Such elastic properties are usually due to hydrogen bonding between residues in  $\beta$  strands of the protein molecules.*

#### **Using Main menu animation tools**

You can now play the movie of the loaded trajectory back and forth, using the animation tools in Figure 5.7.15.

7. By dragging the slider (Fig. 5.7.15), one navigates through the trajectory. The buttons to the left and to the right from the slider panel allow one to jump to the end of the trajectory or go back to the beginning.
8. For example, create another representation for water in the Graphical Representations window: click on the Create Rep button; in the Selected Atoms field, type `water` and hit enter; for Drawing Method, choose Lines; for Coloring Method, select Name.

*This representation of water shows the water droplet present in the simulation. Using the slider, observe the behavior of the water around the protein. The shape of the water droplet changes throughout the simulation, because water molecules follow the protein as it unfolds, driven by the interactions with the protein surface.*

When playing animations, you can choose between three looping styles: *Once*, *Loop*, and *Rock*. You can also jump to a frame in the trajectory by entering the frame number in the window on the left of the “slider” panel.

### **Smoothing trajectories**

9. For clarity, turn off the water representation by double-clicking on it in the Graphical Representations window.

*As you might have noticed, when we play the animation, the protein movements are not very smooth due to thermal fluctuations (as the simulation is performed under the conditions that mimic a thermal bath). VMD can smooth the animation by averaging over a given number of frames.*

10. In the Graphical Representations window, select your protein representation and click on the Trajectory tab. At the bottom, you should see the Trajectory Smoothing Window Size set to zero. As your animation is playing, increase this setting. Notice that the motion gets smoother and smoother as the size of the smoothing window is increased. Commonly used values for this setting are 1 to 5, depending on how smooth you want your trajectory to be.

### **Displaying multiple frames**

We will now learn how to display many frames of the same trajectory at once.

11. In the Graphical Representations window, highlight your protein representation by clicking on it and press the Create Rep button. This creates an identical representation, but note that smoothing is set to zero. Hide the old protein representation.
12. Highlight the new protein representation and click the Trajectory tab. Above the smoothing control, notice the Draw Multiple Frames control. It is set to now by default, which is simply the current frame. Enter `0:10:99`, which selects every tenth frame from the range 0 to 99.
13. Go back to the Draw style tab, and change the Coloring Method to Timestep. This will draw the beginning of the trajectory in red, the middle in white, and the end in blue.
14. We can also use smoothing to make the large-scale motion of the protein more apparent. Go back to the Trajectory tab, and set the smoothing window to 20. The result should look like Figure 5.7.16.


### **Updating selections**

Now, we will see how to make VMD update the selection each frame.


15. Hide the current representation showing all frames, and display only the water representation by double-clicking on it. Change the text in the Selected Atoms field from `water` to `water and within 3 of protein` and hit enter. This will show all water atoms within 3 Å of the protein.
16. Play the trajectory.

*As you can see, although the displayed water atoms may be near the protein for a little while, they soon wander off, and are still shown despite no longer meeting the selection criteria. The Update Selection Every Frame option in the Trajectory tab of the Graphical Representations window remedies this. If the option box is checked, the selection is updated every frame. See Figure 5.7.17.*

17. Quit VMD.

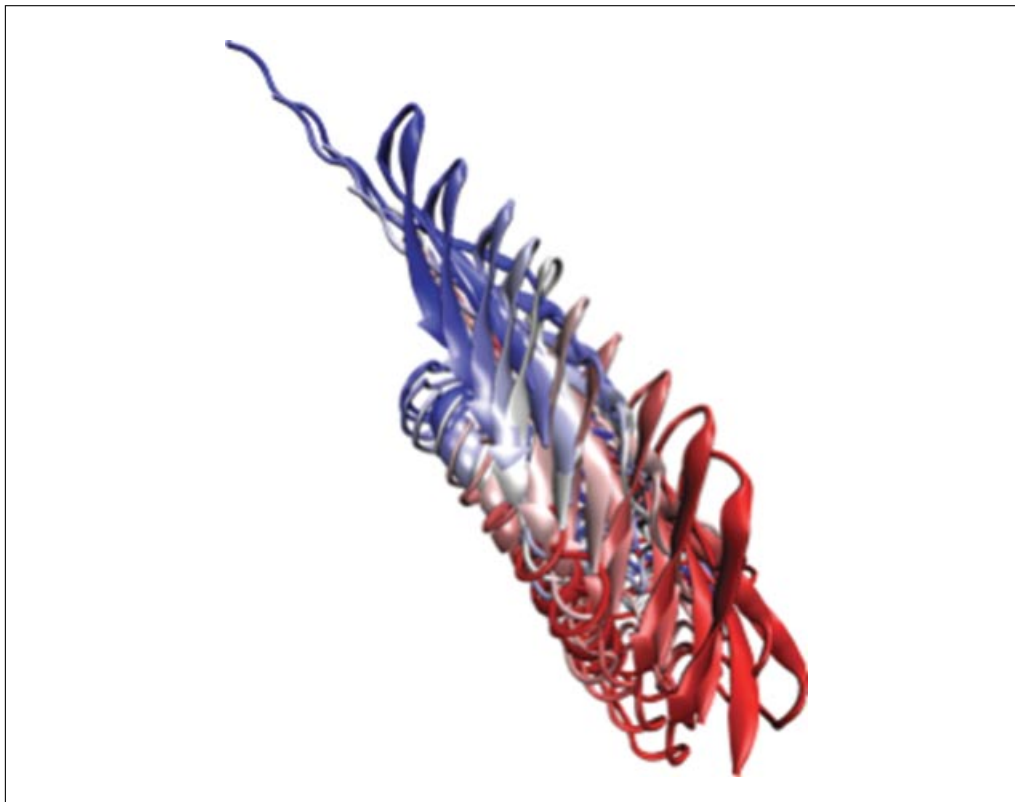

**Figure 5.7.16** Image of every tenth frame shown at once, smoothed with a 20-frame window. For the color version of this figure go to <http://www.currentprotocols.com>.

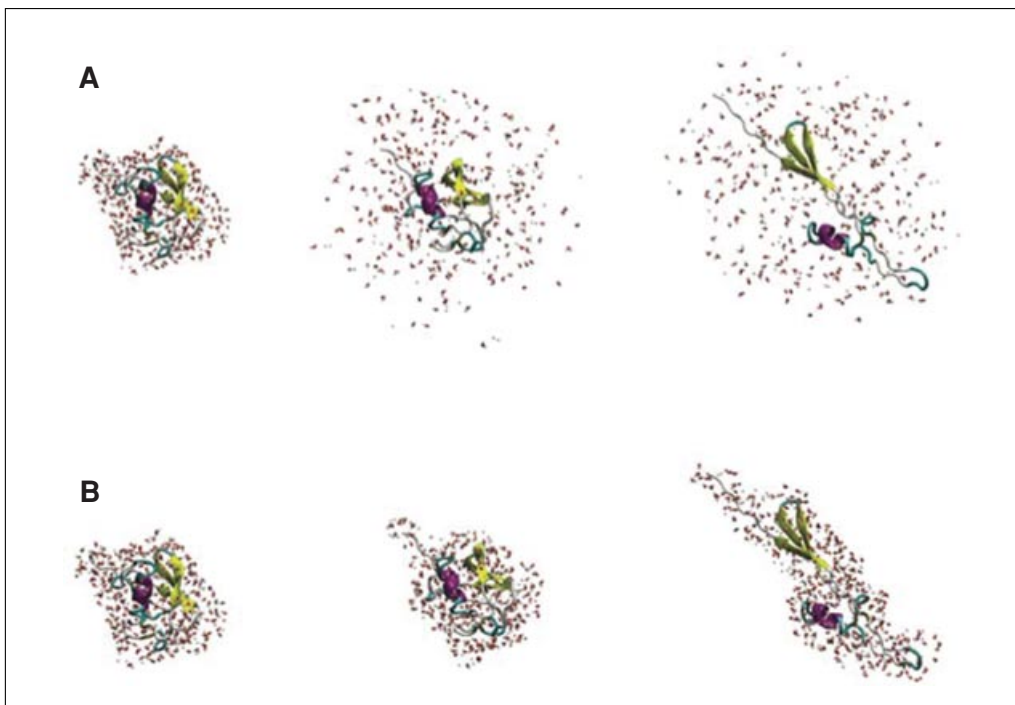

**Figure 5.7.17** Water within 3 Å of the protein, shown for a selection that is not updated (**A**) and for the one that is updated (**B**) each frame. The snapshots shown are (from left to right) for frames 0, 17, and 99. For the color version of this figure go to <http://www.currentprotocols.com>.

## **The Basics of Movie Making in VMD**

The following protocol describes how to make a movie in VMD.

1. Start a new VMD session. Repeat steps 1 to 5 of Basic Protocol 3 to load the ubiquitin trajectory into VMD and display the protein in a secondary structure representation.
2. To make movies, we will use the VMD Movie Maker plugin. In the VMD Main window, go to menu item Extension → Visualization → Movie Maker. The VMD Movie Generator window will appear.

### ***Making single-frame movies***

3. Click on the Movie Settings menu in the VMD Movie Generator window; take a look at the options.

*You can see that in addition to a trajectory movie, Movie Maker can also make a movie by rotating the view point of a single frame. In the Renderer menu, one can choose the type of renderer for making the movie. While renderers other than Snapshot, e.g., Tachyon, generally provide more visually appealing images, they also take longer to render. The rendering time is also affected by the size of the OpenGL window, since it takes more computing time to render a larger image. We will first make a movie of just one frame of the trajectory. Here, we will use the default option, Snapshot. One can also choose the output file format for the movie in menu item Format.*

4. Select the Rock and Roll option in the Movie Settings menu in the VMD Movie Generator window. Set the working directory to any convenient directory of your choice, give your movie a name, and click Make Movie.
5. Once rendering is finished, open and view the movie with your favorite application. This movie setting is good for showing one side of your system primarily.

*If you cannot successfully make movies with VMD, it is possible that you are missing some software required for generating movies. All of the required softwares are freely available, and to find what software you need, please see the VMD Movie Plugin page at <http://www.ks.uiuc.edu/Research/vmd/plugins/vmdmovie/>.*

### ***Making trajectory movies***

6. Now, we will make a movie of the trajectory. In the VMD Movie Generator window, select Movie Settings → Trajectory, give this one a different name, and click Make Movie.

*Note that the length of the movie is automatically set 24 frames per second. For a trajectory, duration of the movie can be decreased, but cannot be increased.*

7. Try out different options in the VMD Movie Generator window. Once you are done, quit VMD.

## **SCRIPTING IN VMD**

VMD provides embedded scripting languages (Python and Tcl) for the purpose of user extensibility. In this section, we will discuss the basic features of the Tcl scripting interface in VMD. You will see that everything you can do in VMD interactively can also be done with Tcl commands and scripts. We will also demonstrate how the extensive list of Tcl text commands can help you investigate molecule properties and perform various types of analysis.

### ***Necessary Resources***

#### ***Hardware***

Computer

## Software

VMD, and a text editor

## Files

1ubq.pdb and beta.tcl, which can be downloaded from <http://www.currentprotocols.com>

## The Basics of Tcl Scripting

Tcl is a rich language that contains many features and commands, in addition to the typical conditional and looping expressions. Tk is an extension to Tcl that permits the writing of graphical user interfaces with windows and buttons, etc. More information and documentations about the Tcl/Tk language can be found at <http://www.tcl.tk/doc>. Let us start with the basic commands.

1. Start a new VMD session. In the VMD Main menu, select Extensions → Tk Console to open the VMD TkConsole window. You can now start entering Tcl/Tk commands here.
2. Try entering the following commands in the VMD TkConsole window. Remember to hit enter after each line and take a look at what you get after each input.

```
set x 10
puts 'the value of x is: $x'
set text 'some text'
puts 'the value of text is: $text'
```

As you can see, the Tcl `set` and `put` commands have the following syntax:

`set variable value` - sets the value of variable

`puts $variable` - prints out the value of variable

Also, `$variable` refers to the value of variable.

3. Try the `expr` command by entering the following lines in the VMD TkConsole window:

```
expr 3 - 8
set x 10
expr -3 * $x
```

The `expr` command performs mathematical operations:

`expr expression` - evaluates a mathematical expression.

4. Entering the following example in the VMD TkConsole window:

```
set result [expr -3 * $x]
puts $result
```

By using brackets, you can embed Tcl commands into others. A bracketed expression will automatically be substituted by the return value of the expression inside the brackets: `[expression]` - represents the result of the expression inside the brackets.

5. Let us calculate the values of  $-3*x$  for integers  $x$  from 0 to 10 and output the results into a file named `myoutput.dat`.

```
set file [open 'myoutput.dat' w]
for {set x 0} {$x <= 10} {incr x} {
puts $file [expr -3 * $x]
}
close $file
```

Here, you have tried the loop feature of Tcl. Tcl provides an iterated loop similar to the “for” loop in C. The `for` command in Tcl requires four arguments: an initialization, a test, an increment, and the block of code to evaluate. The syntax of the `for` command is:

```
for {initialization} {test} {increment} {commands}
```

Take a look at the output file `myoutput.dat`, either by a text editor of your choice, or the command `less` in a terminal window on a Mac or Linux Machine.

### Working with a Molecule Using Tcl Text Commands

Anything that can be done in the VMD graphical interface can also be done with text commands. This allows scripts to be written that can automatically load molecules, create representations, analyze data, make movies, etc. Here, we will go through some simple examples of what can be done using the scripting interface in VMD.

#### Loading molecules with text commands

1. In the VMD TkConsole window, type the command `mol new lubq.pdb` and hit enter.

As you can see, this command performs the same function as described at the beginning of Basic Protocol 1, namely, loading a new molecule with file name `lubq.pdb`.

If you see the error message `Unable to load file 'lubq.pdb' using file type 'pdb'`, you might not be in the correct directory that contains the file `lubq.pdb`. You can use the standard Unix commands in the VMD TkConsole window to navigate to the correct directory.

When you open VMD, by default a vmd console window appears. The vmd console window tells you what’s going on within the VMD session that you are working on. Take a look at the vmd console window. It should tell you a molecule has been loaded, as well as some of its basic properties like number of atoms, bonds, residues etc. The Tcl commands that you enter in the VMD TkConsole window can also be entered in the vmd console window. If you are using a Mac, your vmd console window is the terminal window that shows up when you open VMD.

#### Working with specific parts of a molecule: the `atomselect` command

Many times, you might want to perform operations on only a specific part a molecule. For this purpose, VMD’s `atomselect` command is very useful. The `atomselect` command has the following syntax:

```
atomselect molid “selection command” - creates a new atom selection that includes all atoms described by “selection command”.
```

2. Type `set crystal [atomselect top 'all']` in the Tk Console window.

This command allows you to select a specific part of a molecule. The first argument to `atomselect` is the molecule ID (shown to the very left of the VMD Main window); the second argument is a textual atom selection like what you have been using to describe graphical representations in Basic Protocol 1. The selection returned by `atomselect` is itself a command you will learn to use.

This step creates a selection, `crystal`, that contains all the atoms in the molecule and assigns it to the variable `crystal`. Instead of a molecule ID (which is a number), we have used the shortcut `top` to refer to the top molecule. A top molecule means that it is the target for scripting commands. This concept is particularly important when multiple molecules are loaded at the same time (see Basic Protocol 9 for dealing with multiple molecules in VMD).

The result of `atomselect` is a function. Thus, `$crystal` is now a function that performs actions on the contents of the `'all'` selection.

### **Obtaining and changing molecule properties with text commands**

After you have defined an atom selection, you have many commands that you can use to operate on it. For example, you can use commands to learn about the properties of your atom selection (number of atoms, coordinates, total charge, etc). You can also use commands to change its coordinates and other properties. See VMD User's Guide (<http://www.ks.uiuc.edu/Research/vmd/vmd-1.8.6/ug/>) for an extensive list of commands.

3. Type `$crystal num` in the Tk Console window.

*Passing `num` to an atom selection returns the number of atoms in that selection. Check that this number matches the number of atoms for your molecule displayed in the VMD Main window.*

4. We can also use commands to move our molecule on the screen. You can use these commands to change atom coordinates.

```
$crystal moveby {10 0 0}
$crystal move [transaxis x 40 degree]
```

### **Editing properties of selected atoms**

5. Open the Graphical Representation window by selecting Graphics → Representations... in the VMD Main window. Type in `protein` as the atom selection; change its Coloring Method to Beta and its Drawing Method to VDW. Your molecule should now appear as a mostly red and blue assembly of spheres.

*The "B" field of a PDB file typically stores the "temperature factor" for a crystal structure and is read into VMD's "Beta" field. Since we are not currently interested in this information, we can use this field to store our own numerical values. VMD has a "Beta" coloring method, which colors atoms according to their  $\beta$ -factors. By replacing the Beta values for various atoms, you can control the color in which they are drawn. This is very useful when you want to show a property of the system that you have computed.*

6. Return to the Tk Console window and type `$crystal set beta 0`.


*This resets the "beta" field (which is displayed) to zero for all atoms. As you do this, you should observe that the atoms in your OpenGL window will suddenly change to a uniform color (since they all have the same beta values now).*

*You can obtain and set many atomic properties using atom selections, including segment, chain, residue, atom name, position (x, y and z), charge, mass, occupancy and radius, just to name a few.*

7. In the Tk Console window, type `set sel [atomselect top 'hydrophobic']`.

*This creates a selection, `sel`, that contains all the atoms in the hydrophobic residues.*

8. Let us label all hydrophobic atoms by setting their beta values to 1: type `$sel set beta 1` in the Tk Console window. If the colors in the OpenGL Display do not get updated, go to the Graphical Representations window and click on the Apply button at the bottom.

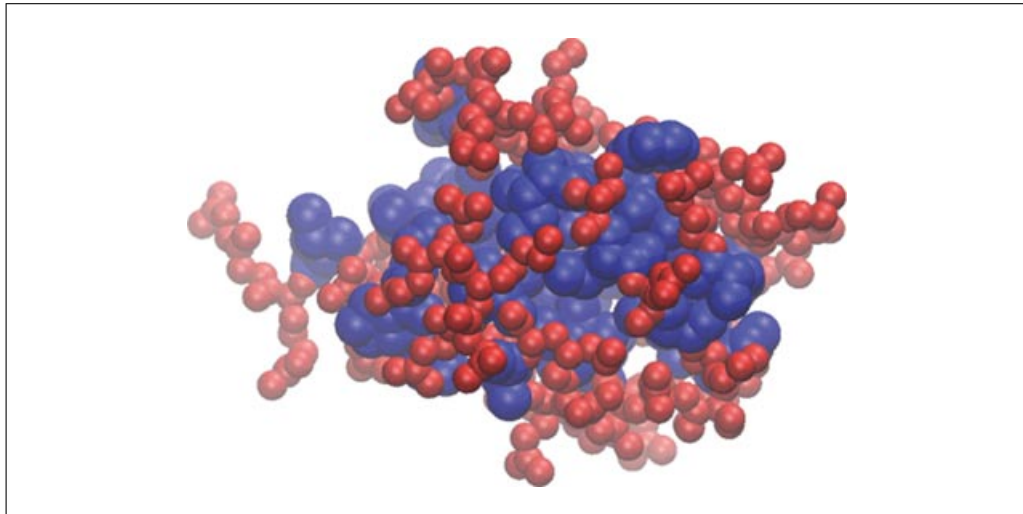

**Figure 5.7.18** Ubiquitin in the VDW representation, colored according to the hydrophobicity of its residues. For the color version of this figure go to <http://www.currentprotocols.com>.

9. You will now change a physical property of the atoms to further illustrate the distribution of hydrophobic residues. In the Tk Console window type `$crystal set radius 1.0` to make all the atoms smaller and easier to see through, and then `$sel set radius 1.5` to make atoms in the hydrophobic residues larger. The radius field affects the way that some representations (e.g., VDW, CPK) are drawn.

*You have now created a visual state that clearly distinguishes which parts of the protein are hydrophobic and which are hydrophilic. If you have followed the instructions correctly, your protein should resemble Figure 5.7.18.*

*Many times in studies of proteins, it is important to identify the locations of the hydrophobic residues, as they often have a functional implication. The method you have just learned is useful in this task. For example, you can easily see that in ubiquitin, the hydrophobic residues are almost exclusively contained in the inner core of the protein. This is a typical feature for small water-soluble proteins. As the protein folds, the hydrophilic residues will have a tendency to stay at the water interface, while the hydrophobic residues are pushed together. This helps the protein achieve proper folding and increases its stability.*

### **The get command**

Atom selections are useful not only for setting atomic data, but also for getting atomic information. For example, if you wish to communicate which residues are hydrophobic, all you need to do is to create a hydrophobic selection and use the `get` command.

10. Try to use the `get` command with your `sel` atom selection to obtain the names of hydrophobic residues:

```
$sel get resname
```

*But there is a problem; each residue contains many atoms, resulting in multiple repeated entries. One way to circumvent this is to pick only the  $\alpha$ -carbons in the selection.*

11. Type the following in the Tk Console window (note, name `CA` =  $\alpha$ -carbons):

```
set sel [atomselect top 'hydrophobic and name CA']
$sel get resname
```

*This should give you the list of hydrophobic residues.*

12. You can also get multiple properties simultaneously. Try the following:

```
$sel get resid
$sel get {resname resid}
$sel get {x y z}
```

*If you want to obtain some of the structural properties, e.g., the geometric center or the size of a selection, the command `measure` can do the job easily.*

13. Let us try using `measure` with the `sel` selection:

```
measure center $sel
measure minmax $sel
```

*The first command above returns the geometric center of atoms in `sel`. And the second command returns two vectors, the first containing the minimum  $x$ ,  $y$ , and  $z$  coordinates of all atoms in `sel`, and the second containing the corresponding maxima.*

Once you are done with a selection, it is always a good idea to delete it to save memory:

```
$sel delete
```

## Sourcing Scripts

When performing a task that requires many lines of commands, instead of typing each line in the Tk Console window, it is usually more convenient to write all the lines into a script file and load it into VMD. This is very easy to do. Just use any text editor to write your script file, and in a VMD session, use the command `source filename` to execute the file. You should have downloaded a simple script, `beta.tcl`, with this unit. We will execute it in VMD as an example. The script `beta.tcl` sets the colors of residues LYS and GLY to a different color from the rest of the protein by assigning them a different beta value, a trick you have already learned in Basic Protocol 6, steps 5 to 9.

In the Tk Console window, type `source beta.tcl` and observe the color change. You should see that the protein is mostly a collection of red spheres, with some residues shown in blue. The blue residues are the LYS and GLY residues in the ubiquitin. Take a quick look at the script `beta.tcl`. Using any text editor of your choice, open the file `beta.tcl`. There are six lines in this file, and each line represents a Tcl command line that you have used before. Close the text editor when you are done.

The `.vmd` file you saved in Basic Protocol 1, step 47, is actually a series of commands. You are encouraged to take a look at that file using a text editor. Hopefully, by the end of this section, you'll understand many of those commands. In fact, you can execute the file in the Tk Console the same way as you execute other script files, i.e., by typing `source myfirststate.vmd` in the Tk Console window.

Many times when you write a script you might want to look up the command for an interactive VMD feature. You can either find it in the VMD User's Guide (<http://www.ks.uiuc.edu/Research/vmd/vmd-1.8.6/ug/>) or conveniently use the console command. Try typing `logfile console` in your Console window. This creates a logfile for all your actions in VMD and writes them in the Console window as command lines. If you execute those command lines, you can repeat the exact same actions you have performed interactively. To turn off logfile, type `logfile off`.

## Drawing Shapes Using VMD Text Commands

VMD offers a way to display user-defined objects built from graphics primitives such as points, lines, cylinders, cones, spheres, triangles, and text. The command that can realize those functions is `graphics`, the syntax of which is `graphics molid command`, where `molid` is a valid molecule ID and `command` is one of the commands shown below. Let us try drawing some shapes with the following examples.

1. Hide all representations in the Graphics Representations window.
2. Let us draw a point. Type the following command in your Tk Console window:

```
graphics top point {0 0 10}
```

*Somewhere in your OpenGL window, there should be a small dot.*

3. Let us draw a line. Type the following command in your Console window (note the “\” in command line means the next line is a continuation of the previous line, hence do not actually type “\” when you enter the following command, and do not start a new line):

```
graphics top line {-10 0 0} {0 0 0} width 5 style\  
solid
```

*This will give you a solid line.*

4. You can also draw a dashed line:

```
graphics top line {10 0 0} {0 0 0} width 5 style\  
dashed
```

*All the objects so far are all drawn in blue. You can change the color of the next graphics object by using the command `graphics top color colorid`. The `colorid` for each color can be found in Graphics → Colors... menu in VMD Main window. For example, the color for orange is “3.”*

5. Type `graphics top color 3` in the Tk Console window and the next object you draw will appear in orange.
6. Try the following commands to draw more shapes:

```
graphics top cylinder {5 0 0} {15 0 10} radius 10\  
resolution 60 filled no  
graphics top cylinder {0 0 0} {-5 0 10} radius 5\  
resolution 60 filled yes  
graphics top cone {40 0 0} {40 0 10} radius 10\  
resolution 60  
graphics top triangle {80 0 0} {85 0 10} {90 0 0}  
graphics top text {40 0 20} 'my drawing objects'
```

7. In your OpenGL window, there are a lot of objects now. To find the list of objects you’ve drawn, use the command `graphics top list`. You’ll get a list of numbers, standing for the ID of each object.
8. The detailed information about each object can be obtained by typing `graphics top info ID`. For example, type `graphics top info 0` to see the information on the point you drew.

9. You can also delete some of the unwanted objects using the command `graphics top delete ID`.

*Using these basic shape-drawing commands, you can create geometrical objects, as well as text, to be displayed in your OpenGL window. When you render an image (as discussed in Basic Protocol 2, steps 19 to 23), these objects will be included in the resulting image file. You can hence use geometric objects and texts to point or label interesting features in your molecule, for example, an arrow (a combination of a cylinder and a cone) can be drawn this way to point at a region of interest of your molecule*

10. Quit VMD.

## WORKING WITH MULTIPLE MOLECULES

In this section, you will learn to work with multiple molecules within one VMD session. We will use the water transporting channel protein, aquaporin, as an example.

### *Necessary Resources*

#### *Hardware*

Computer

#### *Software*

VMD

#### *Files*

1fgy.pdb and 1rc2.pdb, which can be downloaded at <http://www.currentprotocols.com>

### **Molecule List Browser**

Aquaporins are membrane channel proteins found in a wide range of species, from bacteria to plants to human. They facilitate water transport across the cell membrane, and play an important role in the control of cell volume and transcellular water traffic. Many aquaporin protein structures are available in the Protein Data Bank, including a human aquaporin (PDB code 1FQY; Murata et al., 2000) and an *E. coli* aquaporin (PDB code 1RC2; Savage et al., 2003). To practice dealing with multiple proteins in VMD, let us load both aquaporin structures.

### *Loading multiple molecules*

1. Start a new VMD session. In the VMD Main window, choose File → New Molecule. . . The Molecule File Browser window should appear on your screen.
2. Use the Browse. . . button to find the file 1fgy.pdb. When you select the file, you will be back in the Molecule File Browser window. Press the Load button to load the molecule. The coordinate file of human aquaporin AQP1 should now be loaded and can be seen in the OpenGL window.
3. In the Molecule File Browser, make sure you choose New Molecule in the Load files for: pull-down menu on the top. Use the Browse. . . button to find the file 1rc2.pdb and press Load. Close the Molecule File Browser window.


*You have just loaded two molecules. Any number of molecules can be loaded and displayed in VMD simultaneously by repeating the previous step. VMD can load as many molecules as the memory of your computer allows.*

*Take a look at your VMD Main window, which should look like Figure 5.7.19. Within the VMD Main menu you can find the Molecule List Browser (circled in Fig. 5.7.19), which shows the global status of the loaded molecules. The Molecule List Browser displays*

## **BASIC PROTOCOL 9**

### **Modeling Structure from Sequence**

#### **5.7.29**

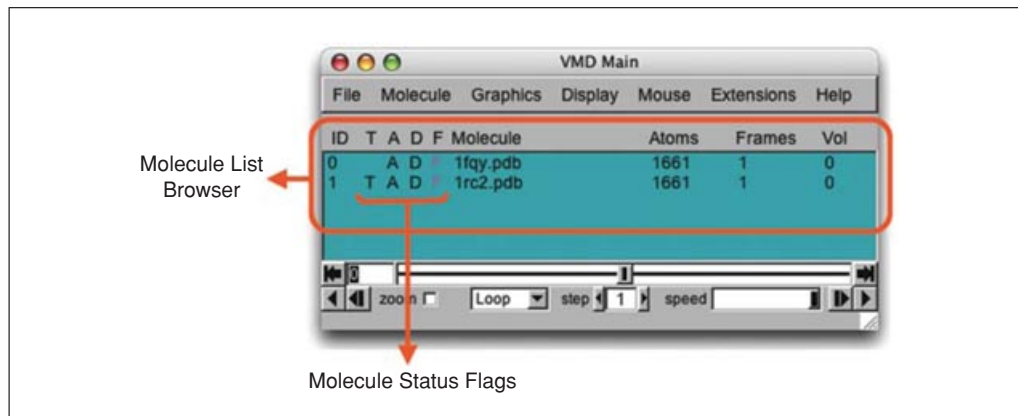

**Figure 5.7.19** The Molecule List Browser.

information about each molecule, including Molecule ID (ID), the four Molecule Status Flags (T, A, D, and F, which stand for Top, Active, Drawn, and Fixed), name of the molecule (Molecule), number of atoms in the molecule (Atoms), number of frames loaded in the molecule (Frames), and the volumetric data loaded (Vol). Let us first start with the Molecule column. By default, the Molecule column displays file names of the molecules loaded in VMD, but you can change the molecule names to recognize them more easily.

### Changing molecule names

- In the VMD Main menu, double-click on `1fqy.pdb` in the Molecule column. A window will pop up with the message Enter a new name for molecule 0:. Type in `human aquaporin`, and click OK (or press enter). In the VMD Main menu, the first molecule now has the name `human aquaporin`.
- Repeat the previous step for the *E. coli* aquaporin by double-clicking the `1rc2.pdb` molecule name, and changing it to `E. coli aquaporin` in the pop-up window.

### Drawing different representations for different molecules

Before we continue exploring other features in the Molecule List Browser, take a look at your OpenGL Display window. You have two aquaporin structures, but since they are both shown in the same default representation, it is difficult to distinguish them. To tell them apart, you can assign them different representations.

- Open the Graphical Representations window via `Graphics → Representations...` from the VMD Main menu. Make sure “0:human aquaporin” is selected in the Selected Molecule pull-down menu on top. Select `NewCartoon` for Drawing Method, and `ColorID → 1 red` for Coloring Method.
- In the Graphical Representations window, select “1:E. coli aquaporin” in the Selected Molecule pull-down menu on top. Select `NewCartoon` for Drawing Method, and `ColorID → 4 yellow` for Coloring Method. Close the Graphical Representations window.

Now, your OpenGL Display window should show a human aquaporin colored in red and an *E. coli* aquaporin colored in yellow.

### Molecule status flags

In your OpenGL Display window, try moving the aquaporins around with your mouse in different mouse modes (rotating, scaling, and translating). You can see that both aquaporins move together. You can fix any molecule by double-clicking the “F” (fixed) flag in the Molecule List Browser on the left of the molecule name.

8. In the Molecule List Browser, double-click on the “F” flag on the left of “human aquaporin” to fix the human aquaporin molecule. Return to the OpenGL Display window and toggle your mouse around. You can see that only the yellow *E. coli* aquaporin moves. Double-click on the “F” flag for human aquaporin again to release it.

*One thing to notice about the “F” flag is that, although it may seem that one molecule has been moved relative to another when one of the molecules is fixed, the difference is only apparent. The internal coordinates of molecules are not changed by the rotation, translation, and scaling motions. To change the coordinates of atoms in a molecule you need to use the text command interface (discussed in Basic Protocol 6, step 4), or by using the atom move picking modes (by choosing Mouse→Move in the VMD Main menu).*

*Other features in the Molecule List Browser include the Molecule ID (ID), Top (T), Active (A), and Drawn (D). Molecule ID is a number (starting from 0) assigned to each molecule when it is loaded into VMD, and permits VMD to recognize each molecule internally. You also refer to molecules by their Molecule IDs in the text command interface. Top flag (T) indicates the default molecule in VMD operations, for example when resetting the VMD OpenGL view and when playing molecule trajectories. There can be only one top molecule at a time. Active flag (A) indicates if the trajectory of the given molecule is updated when using animation tools described in Basic Protocol 3.*

*Finally, Drawn flag (D) indicates if the given molecule is displayed in the OpenGL window. Let us try out the Top and Drawn flags.*

9. Make sure no molecule is fixed. By default, the last molecule loaded in the VMD is the top molecule, so you can check and see that there is a “T” displayed for the *E. coli* aquaporin in the VMD Main menu.
10. Reset the view by pressing the “=” key on the keyboard while keeping the OpenGL Display window active. Note that the yellow *E. coli* aquaporin is now placed in the center of the OpenGL Display window.
11. Switch the top molecule by double-clicking on the empty “T” flag for the human aquaporin molecule in the VMD Main menu. A “T” should appear for the human aquaporin, while the “T” for *E. coli* disappears. Go to the OpenGL Display window and reset the view again. You can see that this time the red human aquaporin is placed in the center of the OpenGL Display window.
12. In the VMD Main menu, try hiding a molecule by double-clicking on its “D” flag. You can display the molecule again by double-clicking its “D” flag again.


### Aligning Molecules with the `measure fit` Command

When you look at your OpenGL Display window, you can see that the two aquaporins are very similar in structure. But it is difficult to detect their slight structural differences as the two proteins are placed apart. We will now try out a very useful Tcl command `measure fit` to align two molecules.

Open the VMD TkConsole window by choosing Extension → TkConsole from the VMD Main menu, and input the following commands:

```
set sel0 [atomselect 0 all]
set sel1 [atomselect 1 all]
set M [measure fit $sel0 $sel1]
$sel0 move $M
```

`measure fit selection1 selection2` – measures the transformation matrix that best aligns the coordinates of `selection1` with the coordinates of `selection2`.

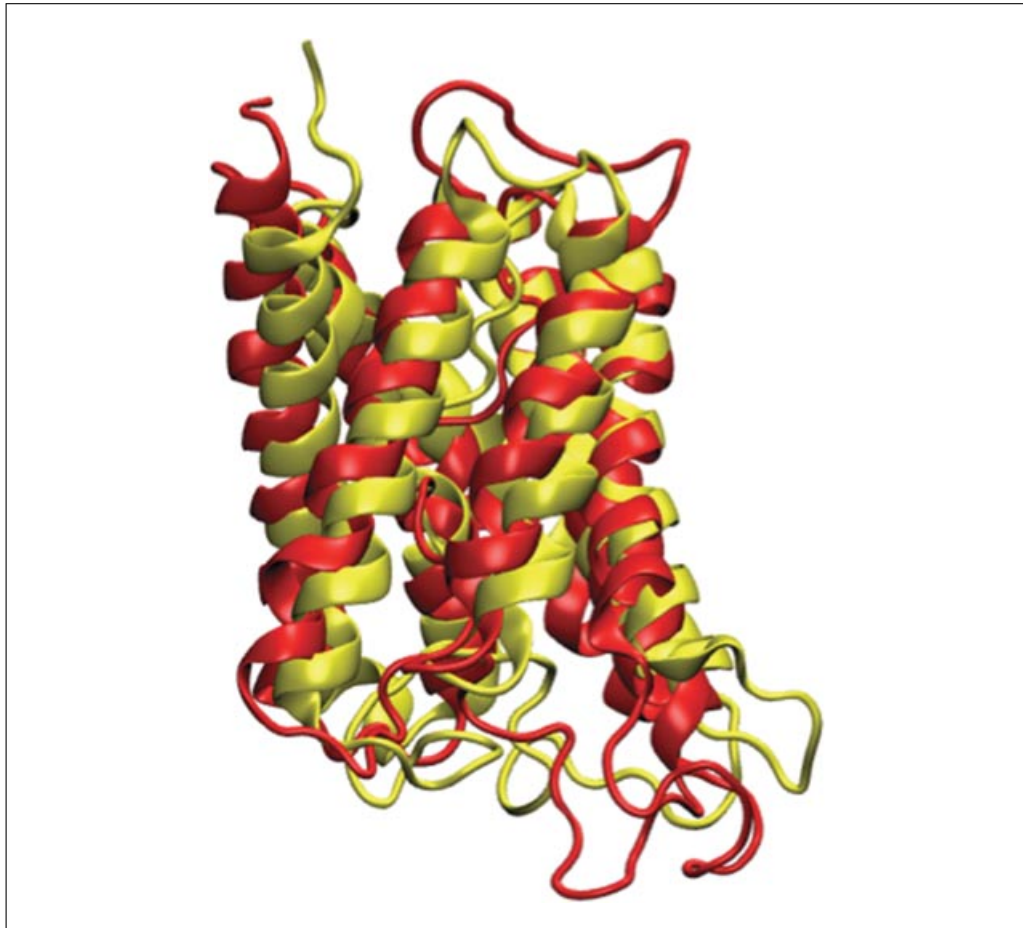

**Figure 5.7.20** Result of the alignment between the two aquaporins using the `measure fit` command. For the color version of this figure go to <http://www.currentprotocols.com>.

As soon as you enter the last command line, you can see that the two aquaporins are now overlapping (Fig. 5.7.20). The  $\alpha$ -helical regions of the aquaporins agree very well, with bigger deviations in the loop regions. Note that the `measure fit` command can only work if two molecules have the same number of atoms. In this case, it is a pure coincidence that the human aquaporin and *E. coli* aquaporin PDB files have the same number of atoms. The `measure fit` command is hence most useful in aligning the same protein in different conformations or different frames of a molecular dynamics simulation trajectory. Generally, to compare the structures of different proteins, one needs to use a different method. A good tool is the VMD MultiSeq plugin, which we will discuss in the following section.

### COMPARING PROTEIN STRUCTURES AND SEQUENCES WITH THE MultiSeq PLUGIN

MultiSeq (Roberts et al., 2006) is a bioinformatics analysis environment developed in the Luthey-Schulten Group at the University of Illinois in Urbana-Champaign. MultiSeq allows users to organize, display, and analyze both sequence and structure data for proteins and nucleic acids, and has been incorporated in VMD as a plugin tool starting with VMD version 1.8.5 (MultiSeq homepage: <http://www.scs.uiuc.edu/~schulten/multiseq>). In this section, you will learn how to compare protein structures and sequences with the VMD MultiSeq plugin. We will again use the water transporting channel protein, aquaporin, as an example.

## Necessary Resources

### Hardware

Computer

### Software

VMD, and a text editor

### Files

1fqy.pdb, 1rc2.pdb, 1lda.pdb, 1j4n.pdb, and spinach\_aqp.fasta, which can be downloaded at <http://www.currentprotocols.com>

## Structure Alignment with MultiSeq

Very often comparing structures of different proteins reveals important information. For example, proteins with similar functions tend to exhibit similar structural features. MultiSeq structure alignment is useful for this reason. We will compare the structures of four aquaporin proteins listed in Table 5.7.8.

### Loading aquaporin structures

1. Start a new VMD session. Open the Molecule File Browser window by choosing the File → New Molecule... menu item in the VMD Main window. In the Molecule File Browser window, use the Browse... button to find and select the file 1fqy.pdb. Press Load to load the molecule.
2. Load the remaining aquaporins, 1rc2, 1lda, and 1j4n. Make sure that each pdb file is loaded into a new molecule. Close the Molecule File Browser window when you have loaded all four molecules. Your VMD Main menu should look like Figure 5.7.21 when all four aquaporins are loaded.

### Aligning the molecules


3. Within the VMD main window, choose the Extension menu and select Analysis → MultiSeq.

*The MultiSeq window (with window name untitled.multiseq showing at the top) should now be open. You may be asked to update some databases in a pop-up window if this is the first time you use MultiSeq. If this is the case, simply click Yes and wait for MultiSeq to finish downloading. When MultiSeq starts, your MultiSeq window should display a list of the four aquaporin protein structures and a list of two nonprotein structures. The nonprotein structures are detergent molecules used in crystallizing the aquaporin proteins, and will not be needed for structure or sequence alignment. You can tell MultiSeq to discard molecules you are not interested in.*

4. In the MultiSeq window, select the 1lda\_X detergent molecule by clicking on it. This will highlight the entire row of 1lda\_X. Remove it from MultiSeq by pressing

**Table 5.7.8** The Four Aquaporins Used in this Section

PDB code	Description	Reference
1fqy	Human AQP1	Murata et al. (2000)
1rc2	<i>E. coli</i> AqpZ	Savage et al. (2003)
1lda	<i>E. coli</i> Glycerol Facilitator (GlpF)	Tajkhorshid et al. (2002)
1j4n	Bovine AQP1	Sui et al. (2001)

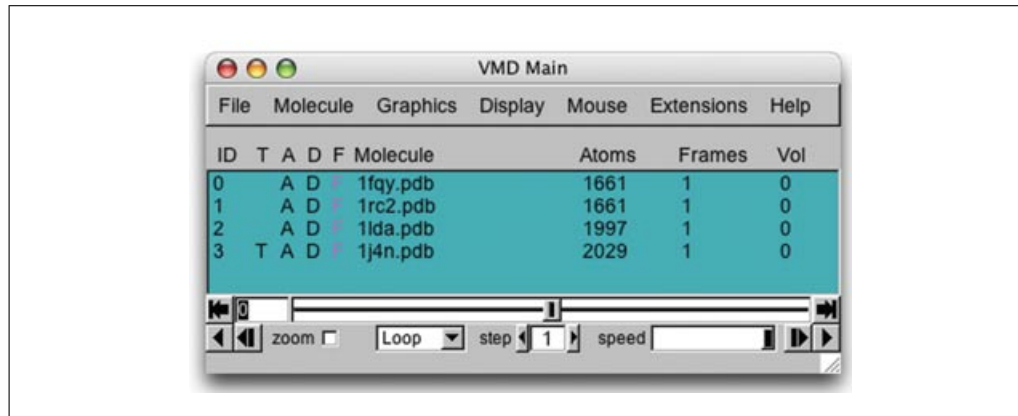

**Figure 5.7.21** VMD Main menu after loading the four aquaporins.

the delete or Backspace key on your keyboard. Do the same to remove the 1j4n\_X detergent molecule.

*MultiSeq uses the program STAMP (Russell and Barton, 1992) to align protein molecules. STAMP (Structural Alignment of Multiple Proteins) is a tool for aligning protein sequences based on three-dimensional structures. Its algorithm minimizes the  $C_{\alpha}$  distance between aligned residues of each molecule by applying globally optimal rigid-body rotations and translations. Note that you can only perform alignments on molecules that are structurally similar; if you try to align proteins that have no common structures, STAMP will fail.*

- In the MultiSeq window, select Tool → Stamp Structural Alignment. This will open the Stamp Alignment Options window.
- In the Stamp Alignment Options window, choose Align the following: All Structures and go to the bottom of the menu and press OK.

*The molecules have been aligned. You can see the alignment both in the OpenGL window and in the MultiSeq window (Fig. 5.7.22). Your alignment in OpenGL window will not immediately resemble Figure 5.7.22. When MultiSeq completes an alignment, it creates a new representation for all the aligned proteins in the NewCartoon representation with the same default coloring method and hides all other representations created previously. Let us give different colors to different aquaporins to distinguish them.*

- Open your Graphical Representations window, and you should see two representations for each molecule, the first one created when VMD loaded the molecule (which is now hidden), and the second one created automatically by MultiSeq. Select “0:1fqy.pdb” in the Selected Molecule pull-down menu on top and highlight the bottom representation by clicking on it. Change the color for this representation by selecting ColorID → 1 red for Coloring Method.
- In the Graphical Representations window, select “1:1rc2.pdb” in the Selected Molecule pull-down menu on top and highlight the bottom representation by clicking on it. Select ColorID → 4 yellow for Coloring Method.
- In the Graphical Representations window, select “2:1lda.pdb” in the Selected Molecule pull-down menu on top and highlight the bottom representation by clicking on it. Select ColorID → 11 purple for Coloring Method.
- In the Graphical Representations window, select “3:1j4n.pdb” in the Selected Molecule pull-down menu on top and highlight the bottom representation by clicking on it. Select ColorID → 12 lime for Coloring Method. Close the Graphical Representations window.


Now your OpenGL window should look similar to Figure 5.7.22, and you can see that the alignment was pretty good as the four aquaporin structures are very similar.

You can get more information about the alignment in the MultiSeq window by highlighting the molecules you wish to compare.


11. In the MultiSeq window, highlight 1fqy by clicking on it.
12. To highlight another molecule without unhighlighting 1fqy, you need to Ctrl-click (or command-click on a Mac) on that molecule. Highlight 1rc2 by clicking on it while holding down the Ctrl key on the keyboard (or the command key on a Mac).

When both 1fqy and 1rc2 are highlighted, you should see at the lower left corner in the MultiSeq window a line of text:  $Q_H:0.6442$ ,  $RMSD:2.3043$ ,  $Percent\ Identity:30.28$ . Note that the values you obtain might be a little different depending on if your MultiSeq database is updated, but they should be close to the ones given here.

The  $Q_H$  value is a metric for structural homology. It is an adaptation of the  $Q$  value that measures structural conservation (Eastwood et al., 2001).  $Q=1$  implies that structures are identical. When  $Q$  has a low score (0.1 to 0.3), structures are not aligned well, i.e., only a small fraction of  $C_\alpha$  atoms superimpose. Along with RMSD and Percent Identity, these numbers tell you that the 1fqy and 1rc2 structures are pretty well aligned. You can repeat the previous step to compare the alignment of other molecules. To unselect a highlighted molecule, Ctrl-click on it again (or command-click on a Mac).

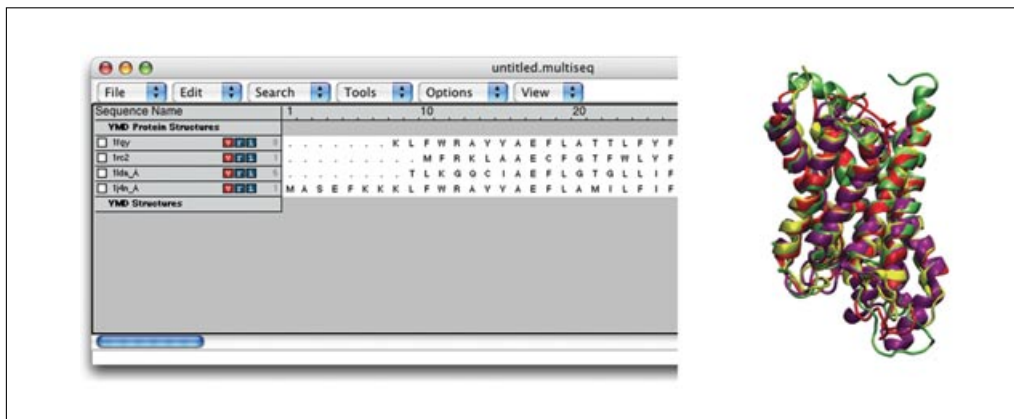

**Figure 5.7.22** The four aquaporins aligned according to their structural similarity. For the color version of this figure go to <http://www.currentprotocols.com>.

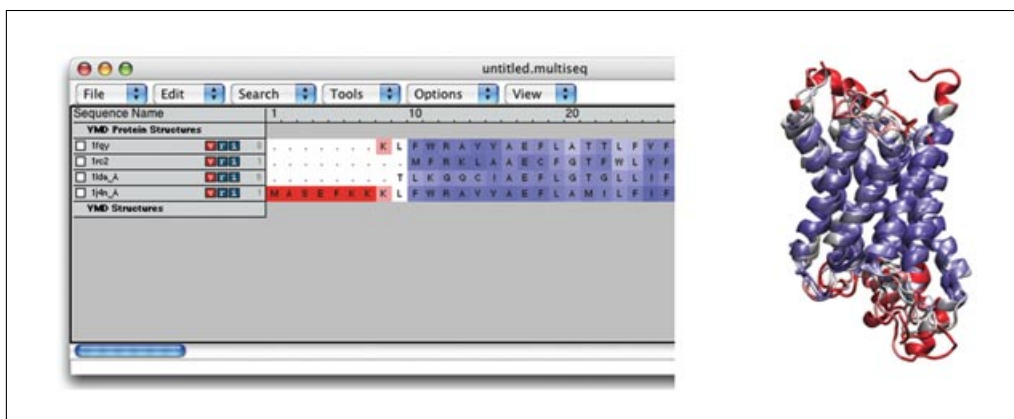

**Figure 5.7.23** Result of a structural alignment of the four aquaporins, colored by  $Q_{res}$ . For the color version of this figure go to <http://www.currentprotocols.com>.

### Coloring molecules according to structural identity

You can also color the molecules according to the value of  $Q$  per residue ( $Q_{res}$ ) obtained in the alignment.  $Q_{res}$  is the contribution from each residue to the overall  $Q$  value of aligned structures.

13. In the MultiSeq window, choose View → Coloring → Qres.

*Look at the OpenGL window to see the impact this selection has made on the coloring of the aligned molecules (Fig. 5.7.23). Blue areas indicate that the molecules are structurally conserved at those points, red areas indicates that there is no correspondence in structure at those points. As you can see, the  $\alpha$ -helices that form the pore are well conserved structurally among the four aquaporins, while there are more structural differences in the less functionally relevant loops.*

### Sequence Alignment with MultiSeq

Besides revealing structural similarities, MultiSeq also allows comparison of proteins based on their sequences. Sequence alignment is often used to identify conserved residues among similar proteins, as such residues are likely of functional importance.

#### Aligning and coloring molecules by degree of conservation


1. In the MultiSeq window, select Tools → ClustalW Sequence Alignment.
2. In the ClustalW Alignment Options window, make sure the Align All Sequences option is checked, and go to the bottom of the window and select OK. Now the four aquaporins have been aligned according to their sequence using the ClustalW tool (Thompson et al., 1994).
3. Let us color the aligned molecules by their sequence similarity. In the MultiSeq window, choose View → Coloring → Sequence identity.

*Now, each amino acid is colored according to the degree of conservation within the alignment: blue means highly conserved, red means low or no conservation. Your MultiSeq window and OpenGL window should resemble Figure 5.7.24.*


*You have now aligned the four aquaporins according to their sequence and identified the conserved residues, found mainly inside the pore (Fig. 5.7.25). Since aquaporin facilitates water transport across the membrane, these conserved residues are most likely the ones that carry out this function.*

#### Importing FASTA files for sequence alignment

Many times the structure of a protein might not be available, but its sequence is. You can analyze a protein in MultiSeq without its structure by loading its sequence information

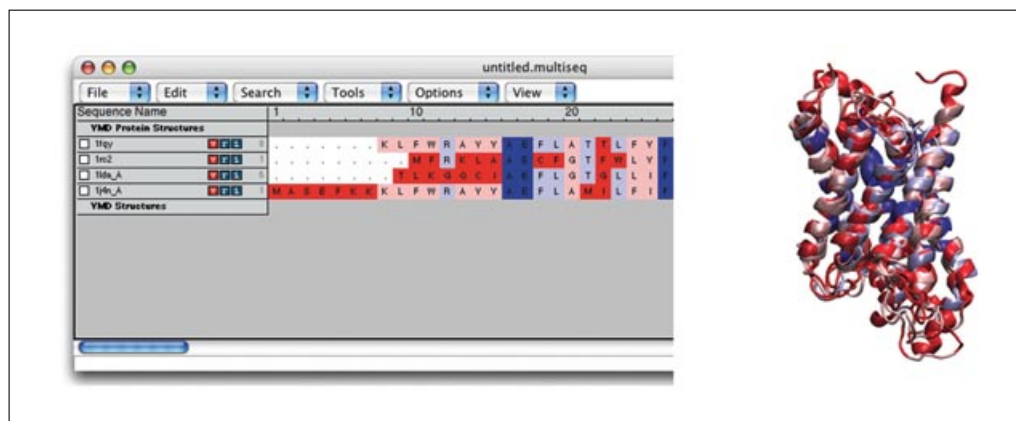

**Figure 5.7.24** Result of a sequence alignment of the four aquaporins, colored by sequence identity. For the color version of this figure go to <http://www.currentprotocols.com>.

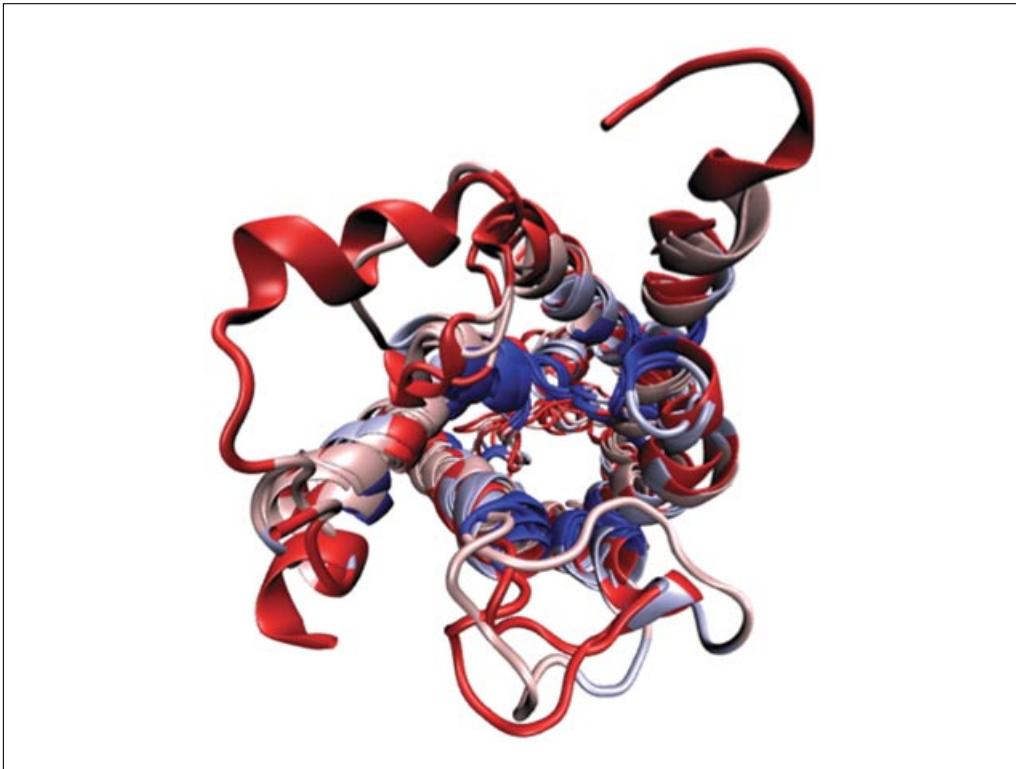

**Figure 5.7.25** Top view of the aligned aquaporins colored by sequence conservation. The conserved residues locate mostly inside the aquaporin pore. For the color version of this figure go to <http://www.currentprotocols.com>.

in the FASTA file format. If you do not have the FASTA file of a protein but you have its sequence, you can create a FASTA file easily with any text editor of your choice.

4. Find the provided FASTA sequence file `spinach_aqp.fasta` and open it with a text editor.

*A FASTA file contains a header that starts with “>” followed by the name of the protein. In the next line is the protein sequence in a one-letter amino acid code. You can create FASTA files similarly in this format. When you create a FASTA file, remember to save it in plain text, and use `.fasta` as the file extension.*

5. Close the text editor when you finish examining `spinach_aqp.fasta`.
6. In the MultiSeq window, select File → Import Data... Select From File in the Import Data window, and press the top Browse button on top to select the file `spinach_aqp.fasta`. Press OK on the bottom of the Import Data window.

*You have now loaded the sequence of a spinach aquaporin into MultiSeq. You can now perform sequence alignment on the spinach aquaporin protein with other loaded aquaporin molecules. Let us try a sequence alignment between a spinach and a human aquaporin.*

7. Click on the checkbox on the left of `spinach_aqp`, and click on the checkbox on the left of `1fqy.pdb`.
8. Open the ClustalW Alignment Options window by selecting Tools → ClustalW Sequence Alignment. Under the Multiple Alignment options on the top, check Align Marked Sequences. Go to the bottom of the window and select OK.

*The sequence of spinach aquaporin is now aligned with the sequence of human aquaporin, and you can check how good the alignment is by obtaining its  $Q_H$  and Sequence Identity values. If you feel that the two molecules are listed too far apart in the MultiSeq window, you can move the molecules by dragging them with your mouse. Also, as you might have noticed, in MultiSeq molecules can be Marked by checking their checkboxes. They can*

also be Selected by highlighting them. You can align only the molecules of your choice by selecting *Align Marked Sequences* or *Align Selected Sequences*, depending if you have marked or highlighted your molecules. This option is available for both structural alignment and sequence alignment.

The structure of spinach aquaporin is actually available (Törnroth-Horsefield et al., 2006), but now that you have learned how to import FASTA sequence data, you can compare the sequences of proteins even if their structures are not resolved yet experimentally.

9. When you finish comparing the sequence of spinach aquaporin with other aquaporins, delete it by clicking on `spinach_aqp` and press delete or Backspace on your keyboard.

## BASIC PROTOCOL 13


### Creating a Phylogenetic Tree with MultiSeq

The Phylogenetic Tree feature in MultiSeq elucidates the structure-based and/or sequence-based relationships between different proteins. Structure-based phylogenetic trees can be constructed according to the RMSD or  $Q$  values between the molecules after alignment; sequence-based phylogenetic trees can be constructed according to the percent identity or ClustalW values (Thompson et al., 1994).

1. Align the structures again by going to the MultiSeq window and selecting Tools → Stamp Structural Alignment.
2. In the Stamp Structural Alignment window, select All Structures, and keep the default values for the rest of the parameters. Press the OK button to align the structures.
3. In the MultiSeq program window, choose Tools → Phylogenetic Tree. The Phylogenetic tree window will open.
4. Select Structural tree using  $Q_H$ , and press the OK button.

*A phylogenetic tree based on the  $Q_H$  values should be calculated and drawn as shown in Figure 5.7.26A. Here you can see the relationship between the four aquaporins, e.g., how the E. coli AqpZ (1rc2c) is related to human AQP1 (1fgy).*

5. You can also construct the phylogenetic tree of the four aquaporins based on their sequence information. Close the Tree Viewer window.
6. You need to perform the sequence alignment again for the four aquaporin proteins. In your MultiSeq window, choose Tools → ClustalW Sequence Alignment, and make sure the Align All Sequences option is checked, and press OK.
7. In the MultiSeq program window, choose Tools → Phylogenetic Tree to open the Phylogenetic tree window again.
8. Select Sequence tree using ClustalW, and press the OK button. A phylogenetic tree based on ClustalW will be calculated and drawn as shown in Figure 5.7.26B.
9. Quit VMD.

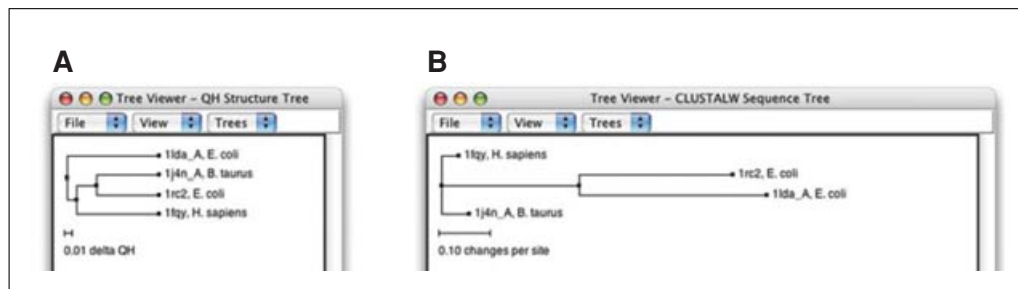

**Figure 5.7.26** (A) A structure-based phylogenetic tree generated by  $Q_H$  values. (B) A sequence-based phylogenetic tree generated by ClustalW.

## DATA ANALYSIS IN VMD

VMD is a powerful tool for analysis of structures and trajectories. Numerous tools for analysis are available under the VMD Main menu item Extension → Analysis. In addition to these built-in tools, VMD users often use custom-written scripts to analyze desired properties of the simulated systems. VMD Tcl scripting capabilities are very extensive, and provide boundless opportunities for analysis. In this section, we will learn how to use built-in VMD features for standard analysis, as well as consider a simple example of scripting.

### *Necessary Resources*

#### *Hardware*

Computer

#### *Software*

VMD, a text editor, and a plotting application

#### *Files*

`ubiquitin.psf`, `pulling.dcd`, `equilibration.dcd` and `distance.tcl`, which can be downloaded at <http://www.currentprotocols.com>

### Adding Labels in VMD

Labels can be placed in VMD to get information on a particular selection, to be used during visualization and quantitative analysis. Labels are selected with the mouse and can be accessed in Graphics → Labels menu. We will cover labels that can be placed on atoms and bonds, although angle and dihedral labelings are also available. In this context, labels for “bonds” or “angles” actually mean distances between two atoms or angles between three atoms, the atoms do not have to be physically connected by bonds in the molecule.

1. Start a new VMD session. Load the ubiquitin trajectory into VMD (using the files `ubiquitin.psf` and `pulling.dcd`). For graphical representation, display protein only, using NewCartoon for drawing method and Structure for coloring method. If you need help, see Basic Protocol 3, steps 1 to 6.
2. Choose the Mouse → Labels → Atoms menu item from the VMD Main menu. The mouse is now set to the mode for displaying atom labels. You can click on any atom on your molecule and a label will be placed for this atom. Clicking again on it will erase the label.
3. We will now try the same for bonds. Choose the Mouse → Label → Bonds menu item from the VMD Main menu. This selects the Display Label for Bond mode.

*We will consider the distance between the  $\alpha$  carbon of Lysine 48 and of the C terminus. In the pulling simulation, the former is kept fixed, and the latter is pulled at a constant force of 500 pN. In reality, polyubiquitin chains can be linked by a connection between the C terminus of one ubiquitin molecule and the Lysine 48 of the next. The simulation then mimics the effect of pulling on the C terminus with this kind of linkage.*

4. Open the TkConsole window by selecting Extensions → Tk Console in the VMD Main menu. We will make a VDW representation for the  $\alpha$  carbons of Lysine 48 and of the C terminus. To find out the indices of these atoms, make a selection including these two atoms by typing in the TkConsole window:

```
set sel [atomselect top ''resid 48 76 and name CA'']
```

**BASIC  
PROTOCOL 14**

**Modeling  
Structure from  
Sequence**

**5.7.39**

5. Get the indices by typing the following line in the TkConsole window:

```
$sel get index
```

This command should give the indices 770 1242.

*Note that the atom numbers of these atoms in the pdb file are 771 and 1243. This is because VMD starts counting atom indices from zero. This is only the case for index, since VMD does not read them from the PDB file. Other keywords, such as residue, are consistent with the PDB file.*


6. In the Graphical Representations window, create a representation for the selection `index 770 1242`, with VDW as drawing method.

7. Now that you can see the two  $\alpha$ -carbons, choose the Mouse  $\rightarrow$  Label  $\rightarrow$  Bonds menu item from the VMD Main menu. Click on each atom one after the other.

*You should get a line connecting the two atoms (Fig. 5.7.27). The number appearing next to the line is the distance between the two atoms in Angstroms. The value of the distance displayed corresponds to the current frame. Try playing the trajectory—you will see that the label is modified automatically as the distance between the atoms changes. Note that the appearance of the line (its color), as well as the appearance of essentially all other objects in VMD, can be changed in Graphics  $\rightarrow$  Colors in the VMD Main menu.*

*The shortcut keys for labels are 1: Atoms and 2: Bonds. You can use these instead of the Mouse menu. Be sure the Open GL Display window is active when using these shortcuts.*

8. The labels can be used not only for displaying, but also for obtaining quantitative information. In VMD Main menu, select Graphics  $\rightarrow$  Labels. On the top left-hand side of the window, there is a pull-down menu where you can choose the type of label (Atoms, Bonds, Angles, and Dihedrals). For now, keep it in Atoms. You can see the list of atoms for which you have made a label.

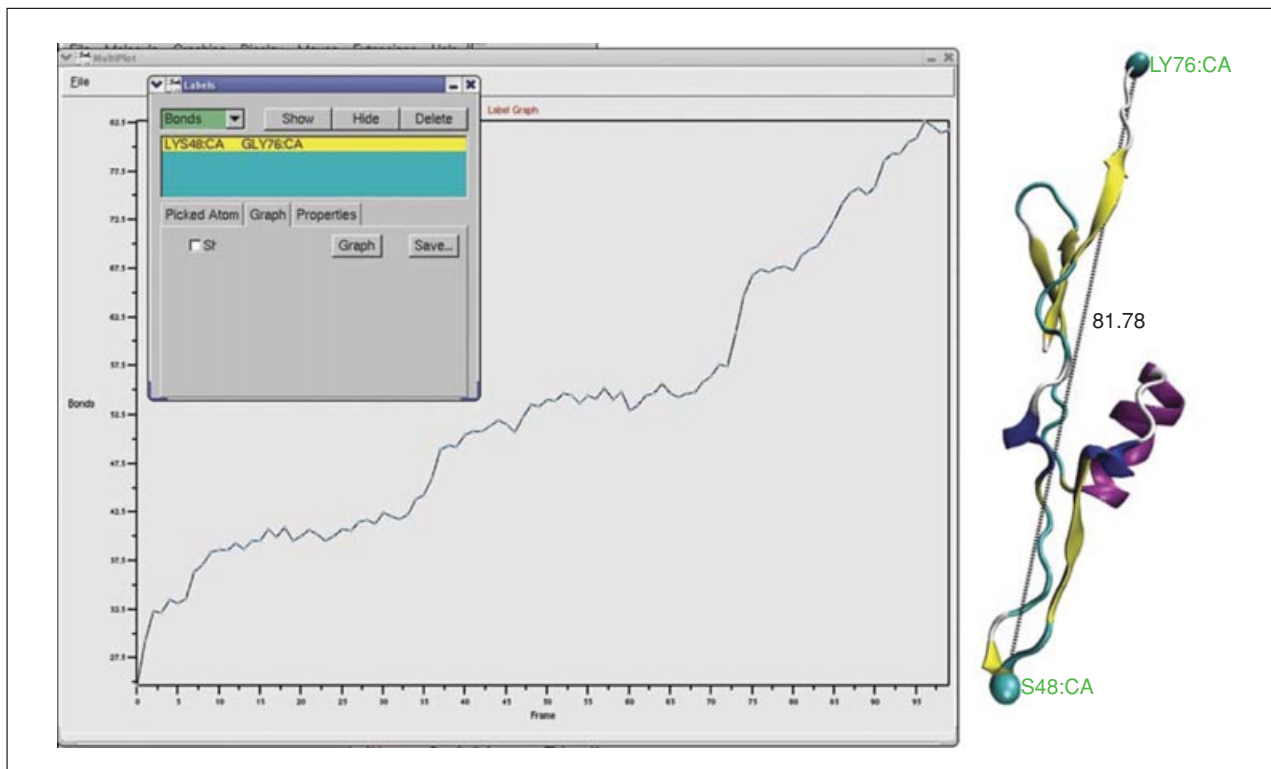

**Figure 5.7.27** Labels in VMD. For the color version of this figure go to <http://www.currentprotocols.com>.

9. Click on one of the atoms. You can see all the information of the atom displayed on the bottom half of the Labels window. This information is useful to make selections; it corresponds to the current frame, and is updated as the frame is changed.
10. You can also delete, hide, or show the atom label by clicking on the corresponding button on the top of the Labels window.
11. In the Labels window, choose label type Bonds, and select the “bond” (distance) you labeled (Fig. 5.7.27). The information given corresponds to only the first atom in the bond, but the number in the Value field corresponds to the length of the bond in Ångstroms.
12. Click on the Graph tab. Select the bond you labeled between atoms 770 and 1242. Click on the Graph button. This will create a plot of the distance between these two atoms over time (Fig. 5.7.27). You can also save this data to a file by clicking on the Save button, and then use an external plotting program to visualize the data.
13. Quit VMD.

### Example of a Built-In Analysis Tool: The RMSD Trajectory Tool

The built-in analysis tools in VMD are available under the menu item Extension → Analysis. These tools each feature a GUI window that allows one to enter parameters and customize the quantities analyzed. In addition, all tools can be invoked in a scripting mode, using the TkConsole window. We will learn how to work with one of the most frequently used tools, the RMSD Trajectory Tool.

In this example, we will analyze RMSD for two trajectories for the same system, `ubiquitin.psf`. One of them is the already familiar pulling trajectory, `pulling.dcd`, and the other is the trajectory of a simulation in which no force was applied to the protein, `equilibration.dcd`.

1. Start a new VMD session. Load the ubiquitin equilibration trajectory into VMD (using the files `ubiquitin.psf` and `equilibration.dcd`).
2. Choose Extension → Analysis → RMSD Trajectory Tool in the VMD Main window (Fig. 5.7.28). The RMSD Trajectory Tool window will show up.


*In the RMSD Trajectory Tool window, you can see many customization options. For the default values, the molecule to be analyzed is `ubiquitin.psf` (the only one loaded). The selection for which RMSD will be computed is all of the protein atoms, excluding hydrogens (since the “noh” checkbox is on). The RMSD will be calculated for each frame with the reference to frame 0. Make sure the Plot checkbox is selected.*

3. Click the Align button.


*This will align each frame of the trajectory with respect to the reference frame (in this case, frame 0) to minimize the RMSD, by applying only rigid-body translations and rotations. This step is not necessary, but is desirable in most cases, because we are interested only in RMSD that arises from the fluctuations of the structure and not from the displacements and rotations of the molecule as a whole. The result of the alignment can be seen in the OpenGL display.*

4. Click the RMSD button in the RMSD Trajectory Tool window. The protein RMSD (in Ångstrom) versus frame number is displayed in a plot (Fig. 5.7.28).

*Over several initial frames,  $RMSD = 0$  because positions of the protein atoms are fixed during that time in the simulation to allow water molecules around the protein to adjust to the protein surface. After that, the protein is released, and the RMSD grows quickly to around 1.5 Å. At that point, the RMSD levels off and remains at ~1.5 Å further on. This is a typical behavior for molecular dynamics simulations. Leveling of the RMSD means*

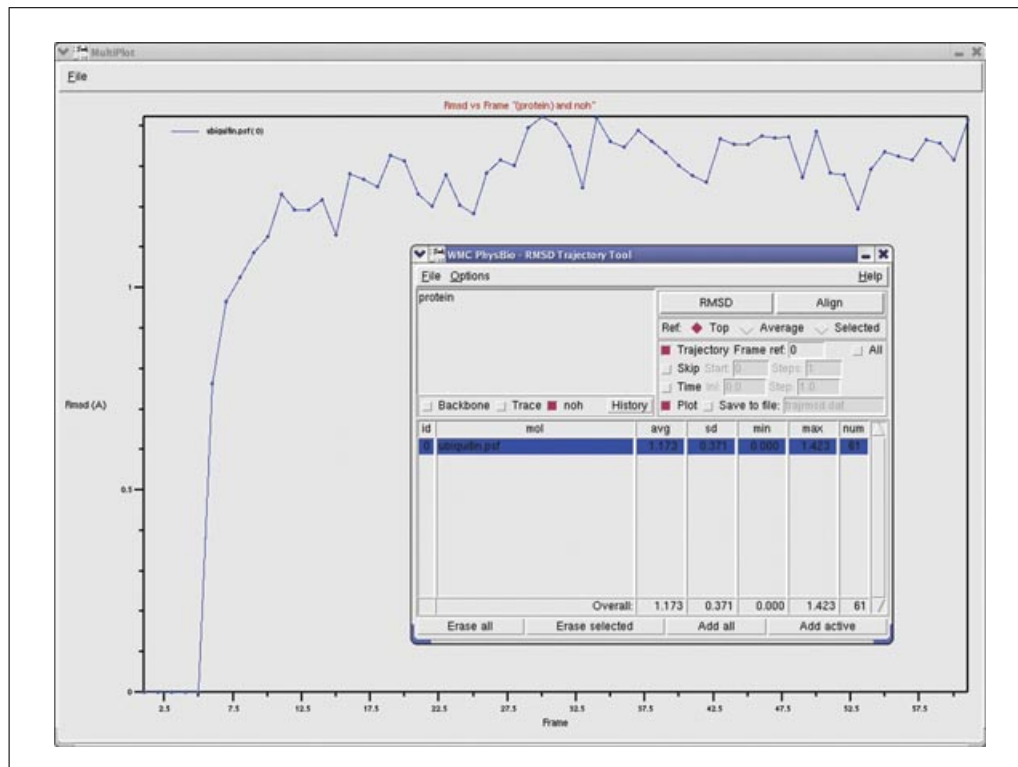

**Figure 5.7.28** RMSD Trajectory Tool. The RMSD is plotted for the equilibration of ubiquitin.

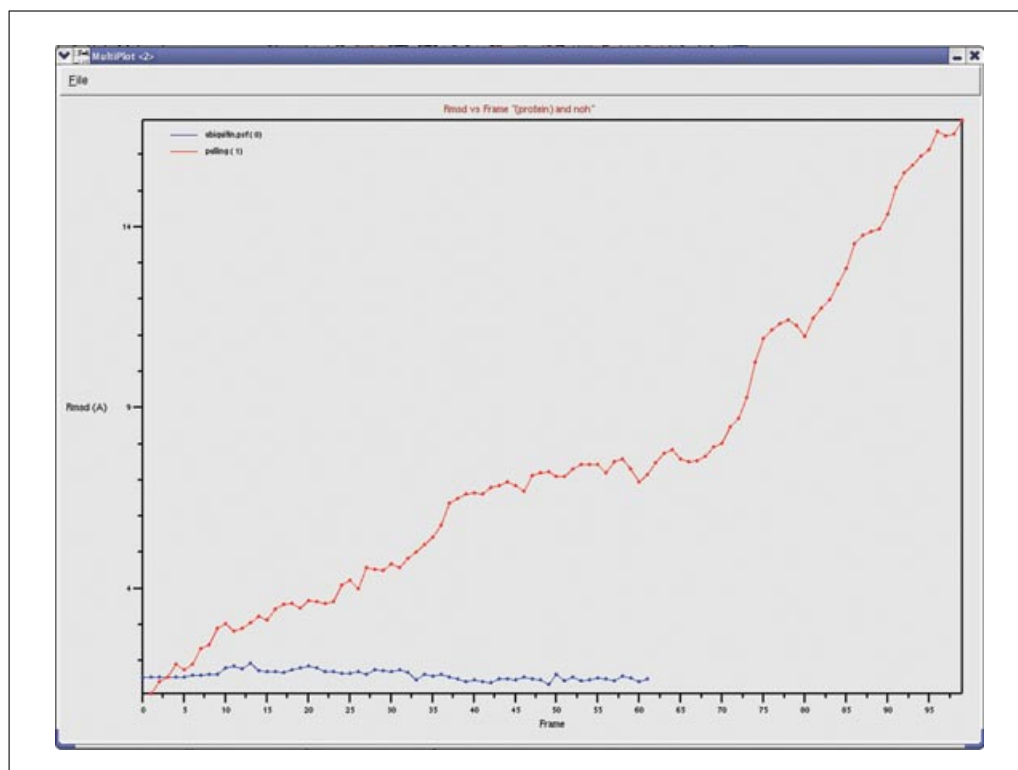

**Figure 5.7.29** RMSD versus time for the equilibration (blue) and pulling (red) trajectories of ubiquitin. For the color version of this figure go to <http://www.currentprotocols.com>.

that the protein has relaxed from its initial crystal structure (which is affected by crystal packing and usually misses some atoms, e.g., hydrogens) to a more stable one. Production molecular dynamics simulations are usually preceded by such equilibration runs, where the protein is allowed to relax; the process is monitored by checking RMSD versus time, and equilibration is assumed to be sufficient when RMSD levels off. The RMSD of 1.5 Å is an acceptable value for most protein simulations. Usually, the deviations from the crystal structure in a simulation are due to the thermal motion and to the relaxation process mentioned; imperfections of the simulation force-fields contribute as well.

5. We will now work with the other trajectory, in which the ubiquitin is pulled apart. Load this trajectory into VMD using the files `ubiquitin.psf` and `pulling.dcd`. Make sure you load `ubiquitin.psf` as a new molecule. You can change the names of the molecules by double-clicking on them in the VMD Main menu (see Basic Protocol 9, steps 4 and 5).
6. In the RMSD Trajectory Tool window, hit the button Add all to update the list of molecules.
7. Click the Align button and then click RMSD button.

*The new graph (Fig. 5.7.29) displays two RMSD plots versus time, one for the equilibration trajectory, and the other for the pulling trajectory. The RMSD for the pulling trajectory does not level off and is much higher than that in the equilibration trajectory, since the protein is stretched in the simulation.*

8. Quit VMD.

### Example of an Analysis Script

In many cases, one requires special types of trajectory analyses that are tailored for certain needs. The Tcl scripting in VMD provides opportunities for such custom tasks. Users commonly write their own scripts to analyze the features of interest. A very extensive library of VMD scripts, contributed by many users, is available online ([http://www.ks.uiuc.edu/Research/vmd/script\\_library/](http://www.ks.uiuc.edu/Research/vmd/script_library/)). Here, we will explore a very simple exemplary script, `distance.tcl`, which computes the distance between two atom selections vs. time and the distribution of the distances.

1. Start a new VMD session. Load the ubiquitin equilibration trajectory (files `ubiquitin.psf` and `equilibration.dcd`).
2. Open the TkConsole window by selecting Extension → Tk Console in the VMD Main menu.
3. In the TkConsole window, load the script into VMD by typing: `source distance.tcl` (make sure that the file `distance.tcl` is in the current folder). This will load the procedure defined in `distance.tcl` into VMD.
4. One can now invoke the procedure by typing `distance` in the TkConsole window. In fact, the correct usage is

```
distance seltext1 seltext2 N_d f_r_out f_d_out
```

where `seltext1` and `seltext2` are the selection texts for the groups of atoms between which the distance is measured, `N_d` is the number of bins for the distribution, and `f_r_out` and `f_d_out` are the file names to where the output distance versus time and distance distribution will be written.

5. Open the script file `distance.tcl` with a text editor. You can see that the script does the following:

- a. Choose atom selections:
 

```
set sel1 [atomselect top '$seltext1']
set sel2 [atomselect top '$seltext2']
```
- b. Get the number of frames in the trajectory and assign this value to the variable `nf`:
 

```
set nf [molinfo top get numframes]
```
- c. Open file specified by the variable `f_r_out`:
 

```
set outfile [open $f_r_out w]
```
- d. Loop over all frames
 

```
for {set i 0} {$i < $nf} {incr i} {
```
- e. Write out the frame number and update the selections to the current frame:
 

```
puts 'frame $i of $nf'
$sel1 frame $i
$sel2 frame $i
```
- f. Find the center of mass for each selection (`com1` and `com2` are position vectors):
 

```
set com1 [measure center $sel1 weight mass]
set com2 [measure center $sel2 weight mass]
```
- g. At each frame `i`, find the distance by subtracting one vector from the other (command `vecsub`) and computing the length of the resulting vector (command `veclegh`), assign that value to an array element `simdata($i.r)`, and print a frame-distance entry to a file:
 

```
set simdata($i.r) [veclength [vecsub $com1 $com2]]
puts $outfile '$i $simdata($i.r)'
```
- h. Close the file:
 

```
close $outfile
```
- i. The second part of the script is for obtaining the distance distribution. It starts from finding the maximum and minimum values of the distance.
 

```
set r_min $simdata(0.r)
set r_max $simdata(0.r)
for {set i 0} {$i < $nf} {incr i} {
set r_tmp $simdata($i.r)
if {$r_tmp < $r_min} {set r_min $r_tmp}
if {$r_tmp > $r_max} {set r_max $r_tmp}
}
```
- j. The step over the range of distances is chosen based on the number of bins `N_d` defined in the beginning and all values for the elements of the distribution array are set to zero.
 

```
set dr [expr ($r_max - $r_min) / ($N_d - 1)]
for {set k 0} {$k < $N_d} {incr k} {
set distribution($k) 0
}
```
- k. The distribution is obtained by adding 1 (`incr ...`) to an array element every time the distance is within the respective bin:
 

```
for {set i 0} {$i < $nf} {incr i} {
set k [expr int(($simdata($i.r) - $r_min) / $dr)]
incr distribution($k)
}
```

1. Write out the file with the distribution:

```
set outfile [open $f.d.out w]
for {set k 0} {$k < $N_d} {incr k} {
  puts $outfile "[expr $r_min + $k* $dr]\
  $distribution($k)"
}
close $outfile
```

6. Now, run the script by typing in the TkConsole window:


```
distance ''protein'' ''protein and resid 76'' 10\
res76-r.dat res76-d.dat
```

*This will compute the distance between the center of the protein and center of the terminal residue 76, and write the distance versus time and its distribution to files res76-r.dat and res76-d.dat.*

7. Repeat the same for the protein's residue 10 by typing in the TkConsole window:

```
distance ''protein'' ''protein and resid 10'' 10\
res10-r.dat res10-d.dat
```

*The data in files produced by the script distance.tcl are in two-column format. Compare the outputs for residue 76 and 10 using your favorite external plotting program (Fig. 5.7.30).*

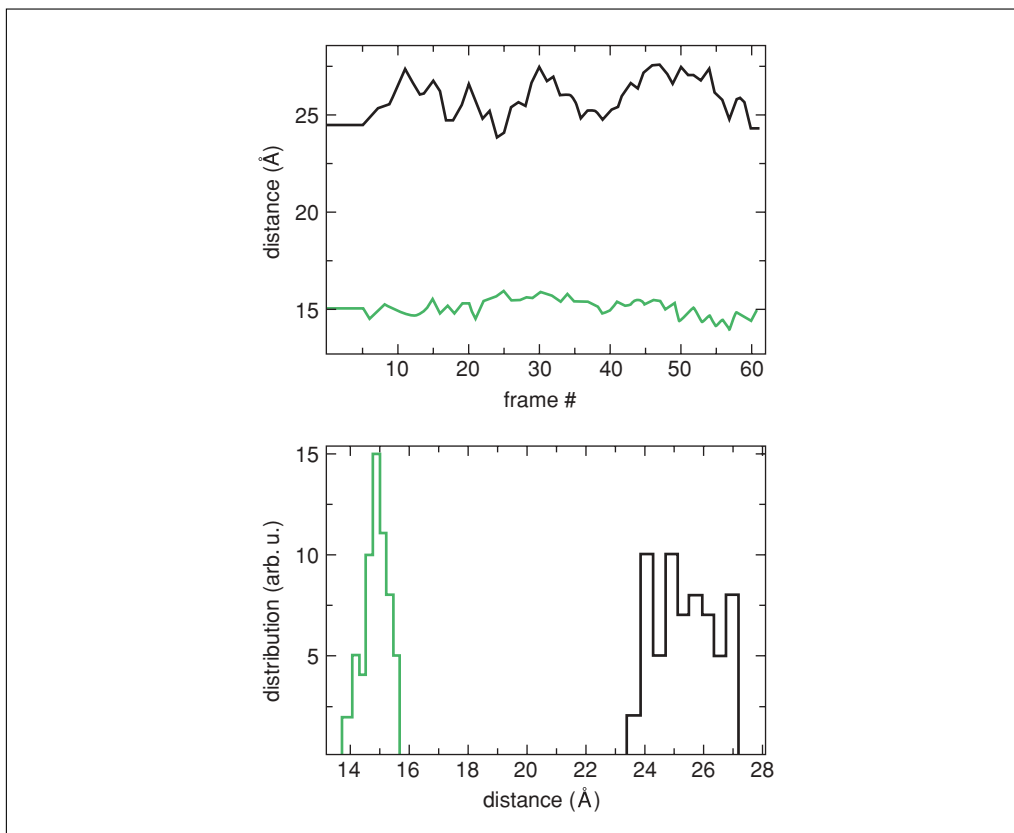

**Figure 5.7.30** Distance between a residue and the center of ubiquitin. The distances analyzed are those for residue 76 (black) and residue 10 (green). For the color version of this figure go to <http://www.currentprotocols.com>.

*Residue 76 is at the protein's C-terminus, which is extended towards the solvent and is quite flexible, while residue 10 is at the surface of the globular part of ubiquitin. The difference in their dynamics with respect to the rest of the protein is immediately obvious when our newly obtained data are plotted (Fig. 5.7.30): the distance of residue 76 from the protein's center is substantially greater than that of residue 10, and the distribution of the distance is noticeably wider due to the flexibility of the C-terminus. This is just a simple example of scripting for the analysis of a trajectory. Similar, but usually much more complex, customized scripts are routinely employed by VMD users to perform many kinds of analysis.*

## 8. Quit VMD.

## COMMENTARY

### Background Information

VMD has been developed by the Theoretical and Computational Biophysics Group at the University of Illinois at Urbana-Champaign. Throughout its development, many features have been added, and user-specific functions can be implemented through embedded scripting languages like Python and Tcl, providing a wide spectrum of tools for the scientific community. Specifically, VMD is most suitable for high-resolution visualization and image rendering, preparation of molecular dynamics simulation systems and analysis of simulation results, and animation of molecular dynamics trajectory. In addition, VMD can also work with volumetric data, and provides a platform for bioinformatics analysis such as protein sequence alignment. What we are able to present in this tutorial only showcases a small part of VMD's capability. But now that you have learned the basics of VMD, you are ready to explore its many other features most suitable for your research. For this purpose, there are many tutorials available that aim at offering a more focused training, either on a specific tool or on a scientific topic. You can find many useful documentations, including the comprehensive VMD User's Guide, at the VMD homepage <http://www.ks.uiuc.edu/Research/vmd/>.

### Critical Parameters and Troubleshooting

Most parameters in VMD can be easily adjusted to suit individual users' needs. For example, when rendering molecules using a representation, as described in Basic Protocol 1, users can adjust the resolution of the representation in the graphical user interface, as well as many other parameters specific to the drawing method of the representation. New users of VMD might find default settings for most parameters are good starting points, but are also encouraged to change the parameters and test the difference. If you have any ques-

tions on using VMD, we encourage you to subscribe to the VMD mailing list [http://www.ks.uiuc.edu/Research/vmd/mailling\\_list/](http://www.ks.uiuc.edu/Research/vmd/mailling_list/).

### Acknowledgments

This tutorial is largely based on the following VMD tutorials, case studies, and user's guides. We hence would like to thank these authors who have provided this tutorial its starting form:

Jordi Cohen, Marcos Sotomayor, and Elizabeth Villa, "VMD Molecular Graphics."

Alek Aksimentiev, John Stone, David Wells, and Marcos Sotomayor, "VMD Images and Movies Tutorial."

Fatemeh Khalili, Elizabeth Villa, Yi Wang, Emad Tajkhorshid, Brijet Dhaliwal, Zan Luthey-Schulten, John Stone, Dan Wright, and John Eargle, "Aquaporins with the VMD MultiSeq Tool."

VMD has been developed by the Theoretical and Computational Biophysics Group at the University of Illinois and the Beckman Institute, and is supported by funds from the National Institutes of Health and the National Science Foundation.

### Citing VMD

The development of VMD is funded by the National Institute of Health. Proper citation is a primary way in which we demonstrate the value of our software to the scientific community, and is essential to continued NIH funding for VMD. The authors request that all published work, that utilizes VMD include the primary VMD citation at a minimum:

Humphrey, W., Dalke, A. and Schulten, K., "VMD - Visual Molecular Dynamics," *J. Molec. Graphics*, 1996, vol. 14, pp. 33-38.

Work that uses softwares or plugins incorporated into VMD should also add the proper citations for those tools. For example, work that uses MultiSeq as introduced in Basic Protocols 11 to 13 should cite:

Roberts, E., Eargle, J., Wright, D. and Luthey-Schulten Z., "MultiSeq: Unifying sequence and structure data for evolutionary analysis," *BMC Bioinformatics*, 2006, 7:382.  
Please see <http://www.ks.uiuc.edu/Research/vmd/allversions/cite.html> for more information on how to cite VMD and its tools.

### Literature Cited

Cruz-Chu, E.R., Aksimentiev, A., and Schulten, K. 2006. Water-silica force field for simulating nanodevices. *J. Phy. Chem. B* 110:21497-21508.

Eastwood, M.P., Hardin, C., Luthey-Schulten, Z., and Wolynes, P.G. 2001. Evaluating protein structure-prediction schemes using energy landscape theory. *IBM J. Res. Dev.* 45:475-497.

Freddolino, P.L., Arkhipov, A.S., Larson, S.B., McPherson, A., and Schulten, K. 2006. Molecular dynamics simulations of the complete satellite tobacco mosaic virus. *Structure* 14:437-449.

Frishman, D. and Argos, P. 1995. Knowledge-based secondary structure assignment. *Proteins* 23:566-579.

Humphrey, W., Dalke, A., and Schulten, K. 1996. VMD—Visual Molecular Dynamics. *J. Mol. Grap.* 14:33-38.

Israelowitz, B., Gao, M., and Schulten, K. 2001. Steered molecular dynamics and mechanical functions of proteins. *Curr. Opin. Struct. Biol.* 11:224-230.

Murata, K., Mitsuoka, K., Hirai, T., Walz, T., Agre, P., Heymann, J.B., Engel, A., and Fujiyoshi, Y. 2000. Structural determinants of water permeation through aquaporin-1. *Nature* 407:599-605.

Phillips, J.C., Braun, R., Wang, W., Gumbart, J., Tajkhorshid, E., Villa, E., Chipot, C., Skeel, R.D., Kale, L., and Schulten, K. 2005. Scalable molecular dynamics with NAMD. *J. Comput. Chem.* 26:1781-1802.

Roberts, E., Eargle, J., Wright, D., and Luthey-Schulten, Z. 2006. MultiSeq: Unifying sequence and structure data for evolutionary analysis. *BMC Bioinformatics*. 7:382.

Russell, R.B. and Barton, G.J. 1992. Multiple protein sequence alignment from tertiary structure comparison: Assignment of global and residue confidence levels. *Proteins* 14:309-323.

Savage, D.F., Egea, P.F., Robles-Colmenares, Y., O'Connell, J.D. III, and Stroud, R.M. 2003. Architecture and selectivity in aquaporins: 2.5 Å X-ray structure of aquaporin Z. *PLoS Biol.* 1:E72.

Sotomayor, M., Vasquez, V., Perozo, E., and Schulten, K. 2007. Ion conduction through MscS as determined by electrophysiology and simulation. *Biophys. J.* 92:886-902.

Sui, H., Han, B.-G., Lee, J.K., Walian, P., and Jap, B.K. 2001. Structural basis of water-specific transport through the AQP1 water channel. *Nature* 414:872-878.

Tajkhorshid, E., Nollert, P., Jensen, M.Ø., Miercke, L.J.W., O'Connell, J., Stroud, R.M., and Schulten, K. 2002. Control of the selectivity of the aquaporin water channel family by global orientational tuning. *Science* 296:525-530.

Thompson, J.D., Higgins, D.G., and Gibson, T.J. 1994. CLUSTAL W: Improving the sensitivity of progressive multiple sequence alignment through sequence weighting, position-specific gap penalties and weight matrix choice. *Nucleic Acids Res.* 22:4673-4680.

Törnroth-Horsefield, S., Wang, Y., Hedfalk, K., Johanson, U., Karlsson, M., Tajkhorshid, E., Neutze, R., and Kjellbom, P. 2006. Structural mechanism of plant aquaporin gating. *Nature* 439:688-694.

Vijay-Kumar, S., Bugg, C.E., and Cook, W.J. 1987. Structure of ubiquitin at 1.8Å resolution. *J. Mol. Biol.* 194:531-544.

Wang, Y., Cohen, J., Boron, W.F., Schulten, K., and Tajkhorshid, E. 2007. Exploring gas permeability of cellular membranes and membrane channels with molecular dynamics. *J. Struct. Biol.* 157:534-544.

Yin, Y., Jensen, M.Ø., Tajkhorshid, E., and Schulten, K. 2006. Sugar binding and protein conformational changes in lactose permease. *Biophys. J.* 91:3972-3985.

Yu, J., Yool, A.J., Schulten, K., and Tajkhorshid, E. 2006. Mechanism of gating and ion conductivity of a possible tetrameric pore in Aquaporin-1. *Structure* 14:1411-1423.

### Supplemental Files

Supplemental files can be downloaded from <http://www.currentprotocols.com> by clicking "Current Protocols" beneath the Bioinformatics head and following the Sample Datasets link.

`1fqy.pdb`

*pdb coordinate file for human aquaporin (Murata et al., 2000)*

`1j4n.pdb`

*pdb coordinate file for bovine aquaporin (Sui et al., 2001)*

`1lda.pdb`

*pdb coordinate file for E. coli GlpF (Tajkhorshid et al., 2002)*

`1rc2.pdb`

*pdb coordinate file for E. coli aquaporin (Savage et al., 2003)*

`1ubq.pdb`

*pdb coordinate file for ubiquitin (Vijay-Kumar et al., 1987)*

`beta.tcl`

*An example tcl script.*

`distance.tcl`

*An example tcl script.*

`equilibration.dcd`

*dcd molecular dynamics trajectory file of an equilibration simulation*

**Modeling  
Structure from  
Sequence**

**5.7.47**

`pulling.dcd`

*dcd molecular dynamics trajectory file of a protein-pulling simulation*

`spinach_aqp.fasta`

*An example fasta protein sequence file.*

`ubiquitin.psf`

*psf structure file for ubiquitin that defines connectivity of atoms*